LANAFORM®

Genial Light

Light Care

Helps to relieve seasonal
affective disorder

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Hept tegen
winterdepressie

GENIAL LIGHT

We would like to thank you for purchasing the GENIAL LIGHT from LANAFORM®. GENIAL LIGHT is an apparatus for dispensing light exposure. It is an efficient and practical means of allowing you to compensate for the lack of natural light. The GENIAL LIGHT provides a very bright light (approximately 10,000 lux) that simulates natural sunshine. You may use the GENIAL LIGHT while continuing other activities such as reading, writing, using the telephone, or even your professional activities.

The additional light provided by the GENIAL LIGHT during the duller seasons increases your exposure to light. It can have a positive influence on your general mood. Having a GENIAL LIGHT means that you can appreciate the benefits of sunlight even on the duller of winter days.

ABOUT LIGHT EXPOSURE

Light and its intensity have a major influence on our vital processes. During the seasons of autumn and winter, solar radiation decreases. Moreover symptoms resulting from a lack of light may appear. Furthermore classic artificial lighting in buildings is not sufficient to compensate for this lack of natural light. The consequences of this situation can lead to seasonal depression, a lack of energy, insomnia or melancholy, as well as the wide range of symptoms covered by the general expression "autumn or winter depression". Light exposure consists of using a light source of high lux value to simulate the natural light of the sun. By undergoing regular exposure to a light source of high lux value over a certain period, the SAD (Seasonal Affective Disorder) symptoms can diminish or even disappear

* A lux is a unit of measurement of light. For example, a sunny summer's day corresponds to approximately 50,000 lux, whereas the typical light intensity within a building is only about 500 lux.

IMPORTANT SAFETY MEASURES


PLEASE READ THESE INSTRUCTIONS BEFORE USING THE APPLIANCE. PLEASE ALSO FOLLOW THESE INSTRUCTIONS CAREFULLY WHILE USING THE APPLIANCE

- Only use this appliance for the purposes described in this manual.
- This appliance is not intended for use by persons (including children) whose physical, sensorial or mental capabilities are impaired, or by persons lacking experience or knowledge, unless they are supervised or instructed in the use of the appliance first by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
- Before connecting the appliance, check that the voltage indicated on the appliance corresponds to the mains supply voltage.
- Please do not use accessories which are not recommended by LANAFORM® or which are not supplied with this appliance.
- It is prohibited to modify this appliance in any way.
- If the power supply lead is damaged, it must be replaced by a similar one available from the supplier or its after-sales service.
- Do not use this appliance if the plug is damaged, if it is not working correctly, if it has been dropped or if it is damaged or if it has fallen into water. If this is the case, have the appliance inspected and repaired by the supplier or its after-sales department.
- Do not transport this appliance by carrying it by its power lead and do not use this lead as a handle.
- Always disconnect the appliance from the mains after use or before cleaning it.
- An electrical appliance should never be left plugged in and unsupervised. Disconnect it from the mains if you are not using it.

- Keep the power lead away from hot surfaces.
- Do not use this product before going to bed. Exposure to light has a stimulating effect which may inhibit sleep.
- Never let anything fall into and never insert anything into one of the openings.
- Do not use this appliance in a humid environment (in a bathroom, near a shower, etc.).
- Never use this appliance in a room where aerosol products (sprays) are used or in a room where oxygen is administered.
- Do not use this appliance under a cover or under a cushion. Excessive heat may cause a fire, electrocution of the individual or other injuries.
- This appliance is intended exclusively for family use. The use of GENIAL LIGHT is therefore not equivalent to a medical treatment.
- Do not use the appliance in a room which is not lit by another source of light, as this can lead to eyestrain. Use this appliance as a supplement for normal lighting (daylight or electric lighting).
- If you feel any pain while using this appliance, stop using it immediately and consult your doctor.
- We recommend that you do not look at the lamp during your session. Pain in the eyes may occur during the first sessions; this will disappear afterwards.
- Consult your doctor before using GENIAL LIGHT if:
 - you are suffering or have suffered from serious depression;
 - you suffer from hypersensitivity to light;
 - you are afflicted with an eye disorder;
 - you are taking certain medications (e.g. certain antidepressants, psychotropic drugs, tablets for malaria) which could have an effect contradicting the use of GENIAL LIGHT;
 - you have high blood pressure.
- Exposure time depends in particular on the distance between you and the GENIAL LIGHT. Refer to the chapter on «Instructions for use» for the recommended distances and exposure times.

COMPONENTS

- 1- 2 x 36W economic PPL fluorescent tubes - colour 84 – giving more than 10,000 lux, with a life span of approximately 8,000 hours
- 2- Reflector hood
- 3- Folding tripod
- 4- Hole for wall-mounting
- 5- Electric power cable
- 6- Cover
- 7- On/Off switch

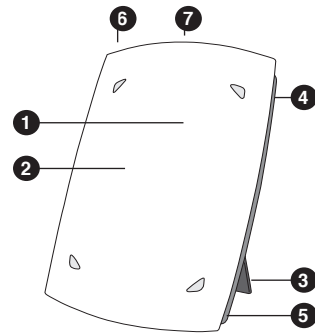INSTRUCTIONS FOR USE

- Place your GENIAL LIGHT on a firm surface, such as a table and adjust the folding tripod. Plug the unit into the power supply and turn on the light by pushing the upper, white button.
- It's recommended to position yourself near the light at a distance of 20 to 50 cm from the lamp depending the time of the session. At a distance of 20 cm, the exposure time should be about 1 hour. At a distance of 50 cm, the exposure time should be about 2 hours.
- During the session, you can read, write, call, etc.

- The exposure time depends upon your position and distance from the light source. The nearer you are to the light, the shorter the session should last.
- The exposure must be regular over a prolonged period of time. During the dullest days of the year, repeat the treatment for 7 consecutive days at least, or more, according to individual requirements.
- After use, switch off the machine by pressing the white button at the top of the unit and unplug it from the power supply.

WARNING

- The GENIAL LIGHT is not a medical instrument.
- Using the GENIAL LIGHT is not a medical treatment.
- We suggest not to look at the lamp during the session. Eye pains or headaches may be experienced during the initial sessions; they will disappear during further sessions.
- If recently you take or have taken medicines which may have a contradictory effect of the use of GENIAL LIGHT, please first consult your doctor;
- If your health deteriorates during your use of GENIAL LIGHT, please stop using the apparatus and consult your doctor.

MAINTENANCE AND STORAGE

Cleaning :

- The lamp must be switched off and unplugged from the power supply before being cleaned.
- Use a damp cloth or sponge to clean the surface of the lamp.
- Never submerge the GENIAL LIGHT in any liquid; and do not use your GENIAL LIGHT in a damp environment.

Storage :

- Put your unit away in a dry and safe place, out of the reach of children.

REPLACING THE LIGHT TUBES

The GENIAL LIGHT operates using two PPL 36W tubes.

- Switch off the lamp and unplug it before replacing the light tubes.
- Remove the screws holding the rear panel of the unit.
- Remove the old tubes and insert two new ones. Replace the screws. N.B.: the old tubes must be disposed of under special conditions. They must not be simply thrown out with domestic waste.

FAQ		
Failure	Possible causes	Solutions
No light	No power	The appliance is not correctly connected or your power supply is defective
	Main switch OFF	Turn on the main switch
One side of the product is out of order	Defective fluorescent bulb	Change the fluorescent bulb
	Fluorescent bulb at end of life	Change the fluorescent bulb
The light continues to flash 20 seconds after startup	Fluorescent bulb incorrectly fitted	Replace the fluorescent bulb

ADVICE REGARDING THE DISPOSAL OF WASTE:

- All the packaging is composed of materials that pose no hazard for the environment and which can be disposed of at your local sorting centre to be used as secondary raw materials. The cardboard may be disposed of in a paper recycling bin. The packaging film must be taken to your local sorting and recycling centre.
- When you have finished using the device, please dispose of it in an environmentally friendly way and in accordance with the law.

LIMITED WARRANTY

LANAFORM® guarantees this product against any material or manufacturing defect for a period of two years from the date of purchase, except in the circumstances described below.

The LANAFORM® guarantee does not cover damage caused as a result of normal wear to this product. In addition, the guarantee covering this LANAFORM® product does not cover damage caused by abusive or inappropriate or incorrect use, accidents, the use of unauthorized accessories, changes made to the product or any other circumstance, of whatever sort, that is outside LANAFORM®'s control.

LANAFORM® may not be held liable for any type of circumstantial, indirect or specific damage.

All implicit guarantees relating to the suitability of the product are limited to a period of two years from the initial date of purchase as long as a copy of the proof of purchase can be supplied.

On receipt, LANAFORM® will repair or replace your appliance at its discretion and will return it to you. The guarantee is only valid through LANAFORM®'s Service Centre®. Any attempt to maintain this product by a person other than LANAFORM®'s Service Centre will render this guarantee void.

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