DE **Haar-Entfernungs-System EPD**

GB **Hair removal system EPD**

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- GB**
- 1 ON/OFF switch
 - 2 Prog + switch
(Program selection forward)
 - 3 + switch (Increase current)
 - 4 Display
 - 5 Belt clip (back page)
 - 6 Eyelets for carrying strap
 - 7 Self-adhesive electrode pads
 - 8 – switch (Reduce current)
 - 9 Prog – switch
(Program selection back)
 - 10 Locking key
 - 11 Battery compartment (back page)
 - 12 Device plug (red/black)
 - 13 Electrode plug
 - 14 Displays the body area
 - 15 Battery display
 - 16 The key lock is activated
 - 17 P , Program number
 - 18 ⌚ , Duration of the application
 - 19 Bar display for current

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**Thank you
very much**

Thank you for your confidence in us and congratulations!

By purchasing the **EPD** hair removal system, you have gained a high-quality product from **MEDISANA**. To derive the most success and to enjoy your **MEDISANA EPD** hair removal system for a very long time, we recommend that you read the following information on the use and care of the device carefully.

IMPORTANT INFORMATION! RETAIN FOR FUTURE USE!

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use.

Should you give this device to another person, it is vital that you also pass on these instructions for use.

You will find the following symbols in the operating manual:

This instruction manual belongs to this device. It contains important information about starting up and operation.

Read the instruction manual thoroughly.

Non-observance of these instructions can result in serious injury or damage to the device.

WARNING

These warning notes must be observed to prevent any injury to the user.

CAUTION

These notes must be observed to prevent any damage to the device.

NOTE

These notes give you useful additional information on the installation or operation.

Safety Information

Tips for your well-being

Do not use the device in the following cases or areas:

- on the front part of the neck (carotid sinus)
- on patients with a pacemaker
- to stimulate the brain region
- in case of pains where the cause is not known
- in case of irritated skin, wounds on the skin or skin infections
- on recent scars
- during pregnancy
- around the eyes or genitals; use on these areas is absolutely prohibited
- in case of metal implants (including the contraceptive coil)
- in areas with reduced sensitivity
- on prosthetic limbs and/or metal dental equipment (for “moustache” treatment)
- on lips; it is absolutely prohibited to position the electrodes on the lips.
- If the patient has any health problems, you may not use the device without prior consultation with the doctor.
- The treatment should not be painful. If you feel pain or feel the application unpleasant, stop it and consult your doctor.
- The skin may not be damaged in any way (skin should not be shaved before the first treatment, in order to prevent any micro-cuts caused by the razor).
- The use of deodorant or perfume is not permitted until 24 hours after epilating.

Always observe the following points

- Only use the device in accordance with its intended use according to the operating instructions.
- The warranty will be invalidated if the device is used for purposes other than those for which it is intended.
- The device may only be used in conjunction with the accessories (electrodes, cable) recommended by the manufacturer.
- Do not use the device if it is not working correctly, if it has been dropped or has fallen into water or has been damaged.
- Do not operate the device in the vicinity of high-frequency transmitters, e.g. microwave and short wave transmitters or radio transmitters. They may impair the device's ability to function.
- Make sure that the user is not in contact with metal objects (e.g. chair or seat).
- Protect the device against damp. If, nonetheless, liquid does enter the device, the batteries must be removed immediately and further use avoided. Please contact your service centre in this case.
- Do not use the device whilst taking a bath or a shower.
- The device may not be used whilst driving a car or whilst operating and controlling machinery.

- This device is not designed to be used by persons (including children) with limited physical, sensory or mental abilities, or by persons with insufficient experience and/or knowledge, unless under observation by a person responsible for their safety, or unless instructed in the use of the device.
- Children must be supervised to ensure that they do not play with the device.
- In case of a malfunction, do not try to repair the device yourself. Not only does the warranty expire, but serious danger can result (fire, electric shock, injuries). Only authorised service offices should carry out a repair.

Informations about the electrodes

- For hygiene reasons, ensure that each user uses his or her own electrode pads.
- Only use self-adhesive electrode pads that stick properly.
- Only position the self-adhesive electrode pads on the places described in this instruction manual.
- Do not place the electrode pads on the body if it is wet.
- Only use the unit on clean, dry skin.
- Before coming into contact with water, remove the electrode pads.
- Do not attach the electrode pads at the front of the neck, the area around the heart, or in the genital area.
- You may not attach the electrode pads on fresh scars, injured or swollen skin, close to thromboses or other vessel problems (e.g. varicose veins).
- Avoid contact with pointed or sharp items.
- Make sure that the two self-adhesive electrodes do not come into contact. This could lead to short-circuiting of the device.
- Do not use damaged or defect leads to the electrode; there is danger of electric shock.
- Do not expose the electrode pads to direct sunlight.
- It is normal for the skin to be slightly reddened after application; the redness will disappear within a few hours. Do not use the unit again until the reddening has disappeared.
- Select shorter treatment times if skin irritations occur over a longer period of treatment.
- In case of more severe skin irritation, stop the treatment and consult a physician.
- Remove the remains of adhesives from the skin with soap and water. Do not use alcohol or solvents.

2.1 Items supplied and packaging

Please check first of all that the unit is complete and is not damaged in any way. If in doubt, do not use the appliance. Send it to a service point.

The following parts are included as standard:

- 1 device **MEDISANA** hair removal system **EPD**
- 2 adhesive electrode pads, approx. 50 x 90 mm (large)
- 4 adhesive electrode pads, approx. 30 x 50 mm (small)
- 2 adhesive electrode pads, approx. 14 x 52 mm (for moustaches)
- 1 connecting cable
- 4 batteries, type AAA (LR03), 1.5 V
- 1 carrying strap
- 1 storage case
- 1 Instructions for use

The packaging can be reused or recycled. Please dispose properly of any packaging material no longer required.

If you notice any transport damage during unpacking, please contact your dealer without delay.

WARNING

Please ensure that polybag packing is kept away from reach of children! Risk of suffocation!

2.2 Hair growth cycles

In order to understand how hair can be permanently removed, one first has to understand how it grows.

- 1) Hair comes from the hair follicles. One square centimetre of human skin can contain between 60 and 500 follicles. A hair can grow from every follicle, but the speed of growth varies for different follicles.
- 2) Every hair passes through a growth cycle, which is divided into the growth phase, the transitional phase and the rest phase. During the growth phase (2 - 7 years), the hair grows approx. 0.3 mm every day. During the transition phase (2 - 3 weeks), the hair root is separated from the nutrient supply, and the hair is pushed towards the surface of the skin. During the rest phase (2 - 4 months), the hair sits loosely in the hair channel and then falls out.
- 3) Hair removal is only effective if the growth of the hair in the hair follicle can be curbed or stopped.
- 4) It is often believed that the hairs appearing one week after epilating are the same hairs as those that were removed, but in this case these are hairs that stem from the neighbouring follicles.
- 5) For this reason, it is important to undergo the treatment one or two times each week, as the first results only appear after a few weeks of treatment and the effect increases through continuous application of the produce.

2.3
How does the EPD hair removal system work?


Normal epilators grip the hairs with their rotating disks, which work like small tweezers, and rip them out along with the roots. This effect lasts for up to four weeks.

The permanent hair removal performed by the **MEDISANA** hair removal system **EPD** goes one step further. A gentle electric impulse from the electrodes attacks the growth cells in the hair roots. The hair root is atrophied. After the hairs have been pulled out, you are free from irritating hairs for a long time.


The **MEDISANA** device is supplied along with special self-adhesive electrode pads for different areas of the body. You can use this pain-free method where you like.

2.4
Electrolysis


During this electrolysis technique, a barely noticeable current is applied to the basis of the hair follicle. A chemical reaction causes the nutrition supply to be cut off. Renewed growth is strongly delayed - so the technology treats not only visible hairs, but also destroys the root bit by bit for long-lasting epilation. It is important that the hair is in its growth phase.

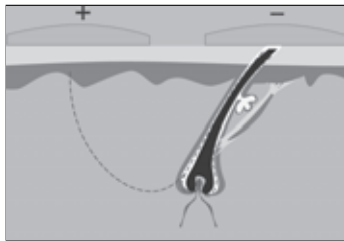Phase 1
Electron flow through the hair

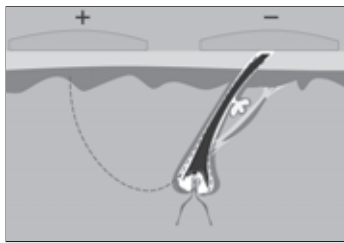Phase 2
Prevention of nutrient supply to the hair follicle

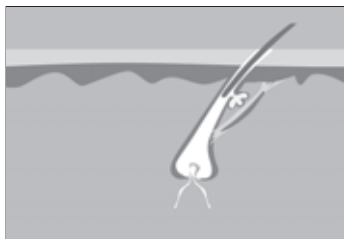Phase 3
The hair that is no longer supplied with nutrients is easy to remove

3.1 Inserting/ changing batteries

Remove the lid from the battery compartment **11** on the back of the device. Insert the four enclosed batteries (type AAA (LR03), 1.5 V). Pay attention to the polarity of the batteries as indicated in the diagram in the battery compartment. Replace the lid of the battery compartment.

When the batteries become weaker, the battery symbol **15** on the display **4** lights up. If the device does not go on after pressing the ON/OFF switch **1**, please insert new batteries.

WARNING BATTERY SAFETY INFORMATION

- Do not disassemble batteries!
- Remove discharged batteries from the device immediately!
- Increased risk of leakage, avoid contact with skin, eyes and mucous membranes! If battery acid comes in contact with skin, rinse the affected area with copious amounts of fresh water and seek medical attention immediately!
- If a battery has been swallowed seek medical attention immediately!
- Replace all of the batteries simultaneously!
- Only replace with batteries of the same type, never use different types of batteries together or used batteries with new ones!
- Insert the batteries correctly, observing the polarity!
- Remove the batteries from the device if it is not going to be used for an extended period!
- Keep batteries out of children's reach!
- Do not attempt to recharge these batteries! There is a danger of explosion!
- Do not short circuit! There is a danger of explosion!
- Do not throw into a fire! There is a danger of explosion!
- Do not throw used batteries into the household refuse; put them in a hazardous waste container or take them to a battery collection point, at the shop where they were purchased!

3.2 Connecting the cable and electrodes

Cable connection

The device is switched off. To connect the device, push the device plug of the cable **12** into the corresponding socket on the device. Plug the red plug (+) into the red socket, and the black plug (-) in the black socket on the device.

Attaching the electrodes

The device is switched off. Remove the carrier film with the two electrode pads from the transport bag. Plug the two electrode plugs at the other end of the cable **13** into the sockets on the electrodes. Leave the electrode pads on the carrier film for protection while doing this.

3.3 Preparing the skin

Before the first treatment, we recommend you clean the skin with an alcohol-free cleanser to remove any remains of grease, make-up etc. Dry the skin thoroughly.

Treatment may only be applied to healthy skin. The treatment may be performed on epilated and non-epilated skin. Hair need not be visible for the electrolysis epilation to be performed (the growth phase begins when the hair is still 4 mm underneath the skin).

3.4 Placing the electrodes

Before placing the self-adhesive electrode pads **7**, please ensure that the device is turned off.

Remove the protective film before attaching the self-adhesive electrode pads. Keep the protective film carefully so that you can stick the electrode pads back onto it after use.

Stick the electrode pads according to the illustration (see chapter **3.5**) to the part of your body where treatment is to be given. Position the electrode pads on the areas of skin where you want to remove the hairs. Press the self-adhesive electrode pads firmly onto the clean and dry skin. If the electrode pads are not attached correctly, no current flow can be established, and the device would switch itself off automatically a few seconds after a program is started (see chapter **3.8**).

WARNING

- **Before you remove the electrode pads or place them elsewhere on your body, always make sure that the device is turned off or turn it off using the ON/OFF switch **1**!**
- **Only then remove the self-adhesive electrode pads carefully by pushing your finger under the edge of the electrode pad and lifting the electrode pad gently thereby. Never remove the self-adhesive electrode pads by pulling the cables.**


CAUTION


- **For safe and hygienic storage stick the electrode pads back onto the protective film.**
- **If the electrode pads **7** have to be replaced, please use only original MEDISANA electrodes pads.**

Once you have attached the adhesive electrode pads on the area of skin to be epilated, you should proceed as described from chapter **3.6**. The size of the electrode pad depends on the surface that is to be treated. Special electrode pads are supplied for the area above the lips. If the body area to be treated allows, we recommend that you use the large electrode pad where possible (50 x 90 mm) (exception: face).

3.5
Examples
of usage

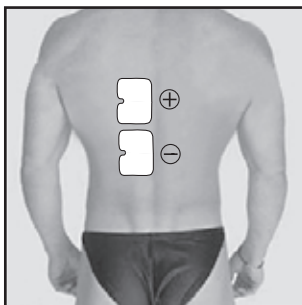
Back

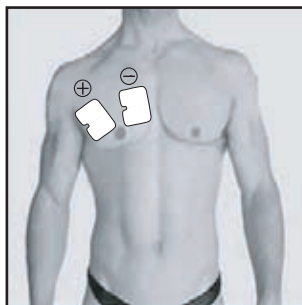Chest


Legs

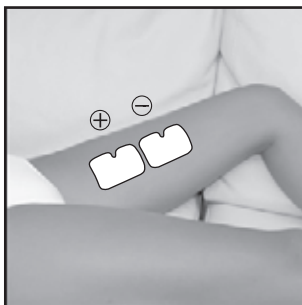
Legs


Bikini zone

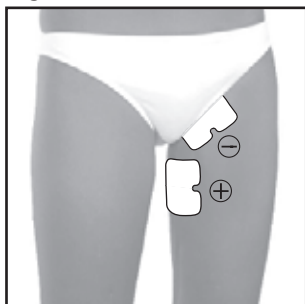Bikini zone

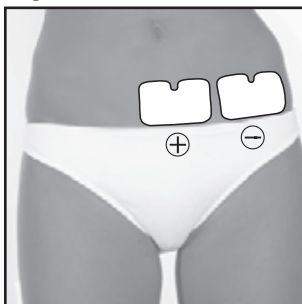
Tummy

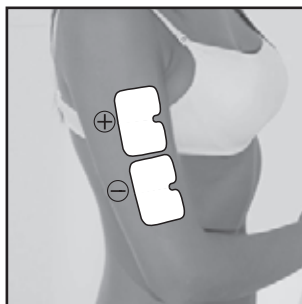Arms

Armpits

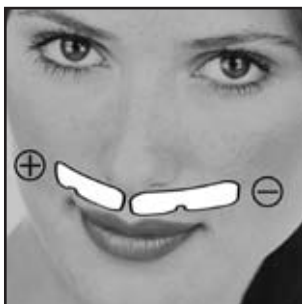Moustache

Face

3.6 Switching on the device

To switch on the device, press the On/off button **1** and hold it for approx. three seconds until you hear an acoustic signal. First all the symbols appear for ca. 1 second in the display **4**. This shows that the device is in working order.

3.7 Programme selection

The EPD device has 12 different programs, 1 - 8 for women and 9 - 12 for men. Use buttons Prog+ **2** and Prog- **9** to select a program. During the selection of the program numbers, the device displays the body areas to be treated **14** and the corresponding program number P **17**. Every press of button **2** or **9** is accompanied by a signal tone.

LIST OF PROGRAMS

No. program	Body area	Duration
<i>WOMAN</i>		
1	Moustache	14 Min.
2	Face	14 Min.
3	Tummy	19 Min.
4	Armpits	14 Min.
5	Bikini zone	12 Min.
6	Legs	14 Min.
7	Arms	12 Min.
8	Chest	10 Min.
<i>MAN</i>		
9	Face	14 Min.
10	Tummy /Back	19 Min.
11	Armpits	14 Min.
12	Legs	14 Min.

3.8 Starting a program

Press the + button **3**. Instead of the program number, the display shows the ⌚ symbol together with the duration of the treatment **18** (10 - 19 minutes) and the bottom bar displaying the strength of current **19**. Unless you regulate the current strength, the program cycle begins with level 1, and the remaining time of the program is displayed in minute increments.

3.9 Increasing/ reducing the current

You can change the strength of the current in 20 steps at any time, either at the beginning of the program or when it is already underway. To increase/reduce the current, press buttons + **3** and - **8**. The bar display only changes when the current is increased by two levels.

Make sure you read the following warnings!

WARNING
INFORMATION AND SAFETY MEASURES ON REGULATING THE CURRENT STRENGTH

The electrolysis effect is created by continuous current. For this reason, you have to select a current that will not irritate or burn the skin.

Always begin with the lowest setting.

First treatment:

We recommend that the first treatment be performed at level 1.

Subsequent treatments:

Set the device to the lowest value at which the patient senses the electricity and heat.

The sensation of heat usually increases during the treatment. Reduce the current strength to make the treatment as pleasant as possible.

As the sensitivity of skin and tissue to continuous current varies for every patient, we recommend that you check the level of redness on the skin after two minutes in the first treatment. To do this, switch off the device and briefly lift the electrodes. The degree of redness is acceptable if it is barely visible. If the redness is easily visible, you have to reduce the current strength. It is normal for the skin to be slightly red at the end of the treatment. This redness goes away after a few hours.

The face and groin ("Bikini zone" program) are especially sensitive and irritable. The current may not be too strong during the first treatments, so that the reaction of the skin can be gauged. In most cases, level 1 should be sufficient.

In most cases, the patient will feel nothing, or only a gentle prickle. The treatment may not become unpleasant or painful. Current and heat must be "barely noticeable".

If even the lowest current level causes the skin to redden seriously, we recommend that the treatment be interrupted and a qualified person be consulted. A slight reddening of the skin is normal.

3.10 Ending the program

The device switches off automatically after the regular program duration has elapsed. If you have to end the program early, switch the device off by pressing the On/Off button **1**.

If the device is switched on and no program is started, the device will switch itself off automatically after 5 minutes.

3.11 Locking the buttons

Once the required current has been set, it is possible to lock the buttons to prevent them from being pressed accidentally. To do this, press and hold the locking button **10** for approx. three seconds. A long signal tone sounds, and the key symbol appears in the display. This means that all buttons apart from the On/Off button are locked. When the buttons are locked, it is not possible to adjust the current or change the program. To unlock the buttons, press and hold the locking button **10** for three seconds again; a long signal tone sounds, the key symbol disappears, and all buttons can be operated again.

3.12 After the treatment

After the treatment, remove the adhesive electrode pads and clean them (see chapter 4.1).

After epilating, you can apply a moisturising lotion to relax the skin.

3.13 Duration of the treatment

The treatment should be performed once a week. The first results appear after a few weeks of treatment, and can be improved through long-term regular application. The duration of the treatment depends on the area being treated, and on the composition of the hairs.

3.14 Frequently asked questions

What has happened if the current returns to 0 immediately after being increased?

The device contains a current control. If the current is increased while the circuit is open, the device resets the current to 0. In this case, you have to ensure that the cables are connected, and that the self-adhesive electrode pads are attached properly to the area to be treated and are not too worn, as this can reduce their conductivity.

Which electrodes may be used?

We recommend using adhesive electrode pads, as they allow better stimulation and are easier to handle. However, the electrode pads have to be replaced when they no longer adhere properly to the skin.

Use only **MEDISANA** original self-adhesive electrode pads.

Where may the electrode pads be positioned?

This instruction manual illustrates the electrode positions for each body zone (see chapter 3.5). Simply follow these instructions.

May this device be used on skin conditions such as psoriasis, hives etc?

No, it is not at all advisable to treat skin with serious skin disease.

4.1 Cleaning and maintenance

- Remove the batteries before cleaning the unit.
- Never use aggressive cleaning agents such as thinning agents, petrol, kerosene, alcohol or chemical products.
- Do not use stiff brushes. Use a soft, damp cloth to clean the device and cables.
- Do not allow any water to get into the unit. Do not use the unit again until it is completely dry.
- Do not expose the device to direct sunlight or to high temperatures and protect it from dust, dirt and damp. Protect the device from jolts, and do not drop it.
- Store the adhesive electrode pads between applications on the protective film in a storage pouch in a cool dry place.
- The self-adhesive electrode pads may not touch or lie above one another. Before every application, make sure that the cables and electrode pads were cleaned properly. If the cables are worn, they may no longer be used; replace them with original parts. Use only self-adhesive electrode pads supplied by the manufacturer.
- The effectiveness of the electrode pads gradually subsides after 20–30 applications due to skin deposits in the gel. To maintain the adhesiveness, hold the electrode pads under cold running water for a few seconds and then shake them dry. Before using the electrode pads again, leave them on the protective film for at least 8 hours.
- Remove the batteries from the unit if you do not plan on using it again anytime in the near future. Otherwise, the batteries may start leaking.
- In case of a malfunction, do not try to repair the device yourself. Not only does the warranty expire, but serious danger can result (fire, electric shock, injuries). Only authorised service offices should carry out a repair.
- The device and accessories are kept in the bag provided.

4.2 Accessories Hair removal system

Ask in your specialist shop or directly at our hotline:

2 self-adhesive electrode pads approx. 50 x 90 mm	Article No.: 87100 9
4 self-adhesive electrode pads approx. 30 x 50 mm	Article No.: 87080 4
2 self-adhesive electrode pads approx. 14 x 52 mm	Article No.: 87110 8
1 connecting cable	Article No.: 85510 8

4.3 Disposal

This product must not be disposed together with the domestic waste.

All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner.

Remove the batteries before disposing of the equipment. Do not dispose of old batteries with your household waste, but dispose of them at a battery collection station at a recycling site or in a shop.

Consult your municipal authority or your dealer for information about disposal.

4.4 Specifications

Name and model	: MEDISANA Hair removal system EPD
Power supply	: 6 V=, 4 x 1,5 V AAA batteries LR03
Output rating	: maximum: 55 V / 5 mA (continuous operation)
Display	: LCD display
Channels	: 1
Impulse form	: continuous
Operating current	: adjustable from 10 mA - 230 mA
Automatic switch off	: after approx. 5 min
Dimensions	: approx. 133 (l) x 83 (w) X 30 (h) mm
Weight	: approx. 140 g
Article No.	: 85500
EAN Code	: 40 15588 85500 9

In accordance with our policy of continual product improvement, we reserve the right for technical and design changes without notice.

Warranty and repair terms

Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

1. The warranty period for **MEDISANA** products is three years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

MEDISANA AG

Itterpark 7-9

40724 Hilden

Germany

Tel.: +49 (0) 2103 / 2007-60

Fax: +49 (0) 2103 / 2007-626

eMail: info@medisana.de

Internet: www.medisana.de

The service centre address is shown on the last page.

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D&L Marketing Oy
ALERE Electronics AB