MEDISANA®

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Foot and back massager FRM

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Thank you very much

Thank you for your confidence in us and congratulations on your purchase! You have purchased a quality product from **MEDISANA**, one of the world's leading companies in the healthcare market. To ensure the best results and long-term satisfaction with your **MEDISANA FRM** massager, we recommend that you read the following operating and maintenance instructions carefully.

1.1 Important

Read the following safety information carefully before using the device and keep these instructions in a safe place for later use.

Note the following:

- Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your mains supply.
- Only use this unit according to the intended purpose and as directed in the instructions for use. The device is not intended for commercial use.
- Never use the unit outdoors.
- Do not sit or stand on the unit; it is not designed to support your full body weight.
- Never cover up the ventilation holes on the bottom of the device.
- Do not put any foreign objects into any of the openings. Avoid contact with pointed or sharp objects. The warranty is null and void if it is not used for the intended purpose.
- If the unit is to be used by or on people who are handicapped in any way or in the presence of children, this must always be done under supervision. Medical equipment is not designed to be used as a toy!
- Do not attempt to repair the unit yourself in the event of any defects. Not only does this render the warranty null and void but you can also run the risk of serious hazards (fire, electric shock, injury). Only have repairs carried out by authorised service points.

Note the following in particular concerning electrical equipment:

- Do not reach for the unit if it has fallen into water. First pull out the mains plug at once.
- Do not use any electrical equipment while in the bath or shower. Do not store or keep any electrical equipment in places where it could fall into the bathtub or washbasin.
- Do not use the unit under the bedclothes or under cushions.
- Keep the mains cable away from hot surfaces.
- Never use the unit in rooms in which aerosols (sprays) are being used or where pure oxygen is being administered.
- Never use the unit if the power cable has been damaged, the mains adapter is defective, if the unit does not work perfectly, if it has fallen into water or has been damaged.
- Never carry, pull or rotate the unit by the mains cable and never pinch the cable between other objects.

1.1 Important

- Do not allow the unit to become overheated due to excessively long use.
- Only pull the plug out of the socket when you have dry hands.

1.2 Important recommendations concerning your health

Every individual reacts differently to a massage. Depending on the individual, some prefer a gentle massage, some find a more intensive massage more pleasant, with or without heat. After trying out the **MEDISANA FRM** massager a few times, you will soon find the ideal setting that achieves the desired relaxation for you. The massager is extremely easy to use. However, you should bear in mind the following tips:

- If you are pregnant, suffer from diabetes or have a pacemaker, contact your doctor before using the unit.
- You should also contact your doctor before initial use of the unit, should you have any reservations about your personal tolerance to treatment or have any health-related questions.
- Do not use the unit on any parts of the body that are swollen, injured or inflamed.
- If you experience pains in the muscles and joints that persist after a long period of time, please inform your doctor. Persistent pains can be symptomatic of serious illness.
- We would advise patients suffering from varicose veins against using the unit.
- Should you experience pain or discomfort during treatment, stop treatment immediately.
- If you experience unexplained pain or swelling in the legs or feet, or suffer from a muscle injury, you should consult your doctor before using the unit. Stop the massage if you experience pain during or swelling after treatment.

2.1 Scope of delivery

Please check first that the unit is complete.

The package contains:

- The FRM massager
- Operating instructions

2.2 The MEDISANA FRM Massager: a wide range of applications

Our feet and back are under great strain each and every day, making them especially vulnerable to tenseness, muscle pains and hardened muscles. The MEDISANA FRM massager lets you give your feet and back a treat, helping them to relax after a long day.

The massager offers you a range of options for achieving this:

- Upper and lower back massage
- Foot massage
- Reflex zone massage of the feet
- Temperature cycle treatment

You have the option to combine each massage with hot/cold treatment. The unit also allows you to gently heat the vibro-node surfaces. This promotes relaxation and has a soothing effect on the muscles.

The unit is also equipped with acupressure massage nodes for a soothing reflex zone massage. This type of massage is based on the knowledge that each organ and part of the body can be has its own specific reflex zone on the soles of the feet. When these zones on the soles of the feet are stimulated, nerve impulses are transmitted to the brain and relayed from there to the organs or other parts of the body. The reflex zones can be used to stimulate or soothe individual organs, stimulate the circulation and relieve pains.

As well as massage, the unit is also suitable for heat treatment, during which hot and cold stimuli transmitted through the skin. This stimulates the tissue beneath the skin, aiding the body's own regenerative capabilities.

3.1 Foot massage and reflex zone massage

Foot massage

Place the **MEDISANA FRM** massager on the floor in front of a comfortable chair and plug it into the mains outlet. Sit down in a comfortable position and rest your feet on the unit.

You can adjust the angle of the unit to your own particular requirements by turning the two adjustable feet on the underside of the unit.

Pressing the massage button **2** once activates the vibro-node massage surface **5** and the acupressure nodes **6** at level **1**. The LED next to "LO" lights up, indicating that low level vibration is active. A further press of the massage button **2** switches to level **2** vibration, and the LED indicator changes from "LO" to "HI". Pressing the massage button **2** one more time switches of the vibration.

You may heat the massage surface by pressing the heat function button **1**. The heat function may be activated on its own or in combination with vibration. The LED next to the button **1** indicates that the heat function is active. You may switch off the heat by pressing the button **1** one more time.


Switch off both functions (vibration and heat) before removing the plug from the mains outlet.

Reflexzone massage

The acupressure nodes **6** are ideal for targeting tenseness in specific areas of the body or just generally improving your feeling of well being. Use the following reflex zone map as a guide.

Concentrate your massage on those areas of the soles of the feet that correspond to the part of the body you wish to treat. Stimulate the required area for three to five minutes.

Reflex zone map:

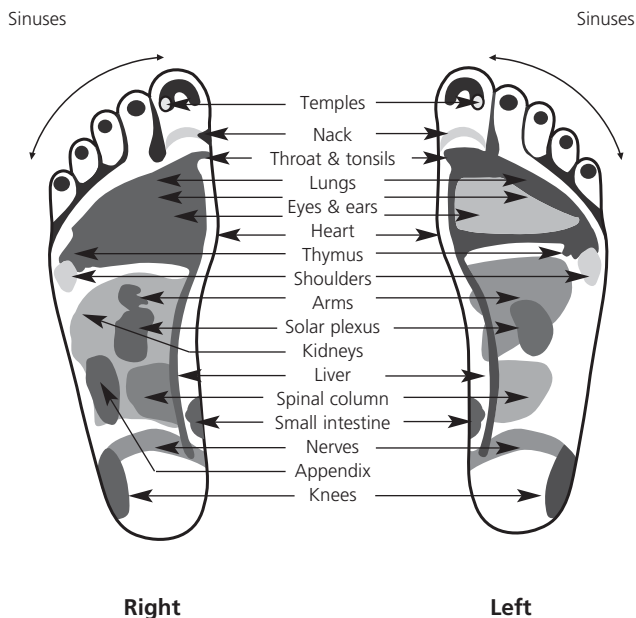3.2 Back massage

Upper back massage

- Turn the adjustable feet anticlockwise until you can feel resistance. You can now hang the unit over the back of a chair. A chair with a high back, such as an office chair, is best suited to this purpose.
- Plug the unit into the mains, making sure that the lead is not stretched. The control buttons ② and ① are operated in the same way as described in the instructions for foot massage.
- You can use this method primarily for treating hardened muscles and soothing tenseness in the neck and upper back region.
- Once you have completed the massage/heat treatment, keep pressing buttons ② and ① until all the LEDs go out and the massager is switched off. Only then should you unplug the unit from the mains outlet.

Lower back massage

- Place the massager vertically on the seat of a chair, with the mains connector ⑧ facing upwards, such that the unit is leaning against and supported by the back of the chair. Plug the unit into the mains, making sure that the lead is not stretched.

- The control buttons are operated in the same way as with a foot massage.
- Sit down in the chair, resting your back against the unit for treating tenseness and hardened muscles in the lower back region.
- Make sure all functions are switched off before unplugging the unit from the mains.

3.3 Temperature cycle treatment

- Place the massager on the floor in front of a comfortable chair and plug it into the mains outlet.
- Before switching on the heated surface, use the hot/cold switch **4** to select the required temperature range. Pushing the hot/cold switch **4** to the right selects heat treatment, while pushing the switch **4** to the left selects cold.
- The raised heat treatment surface will start to heat up or cool down once you press button **3** on the unit. A further press of the button **3** increases the temperature range to level **2**. At the same time, the LED indicator changes from "LO" to "HI". Press the button **3** one more time to switch off the heat treatment. The LED next to the button **3** will go out.
- Alternating hot and cold treatment for three to five cycles stimulates the soles of the feet and promotes better circulation. If you wish, the treatment may be combined with the massage function.
- You should not carry out more than five heat cycle sessions on any one day.

**3.4
Duration of
treatment**

In general, you can allow approximately 10 to 15 minutes per massage session, one to two times a day.
However, you should make sure you allow the unit to cool down completely after each use before using it again!

**4.1
Cleaning and
maintenance**

- Before cleaning the unit, make sure it is switched off and disconnected from the mains. Clean the unit once it has cooled down using a soft, damp cloth.
- Do not use any strong detergents for cleaning, e.g. scouring cream or other cleansing agents, as they may damage the surface of the unit.
- Store the unit in a cool, dry place.
- Take care when rolling up the mains lead to avoid damaging it.

**4.2
Disposal**

Dispose of your old electrical unit in an environmentally acceptable manner when it has reached the end of its service life!
Packaging is re-usable or may be re-cycled; metal parts must be disposed of as scrap metal; plastics, electrical and electronic components must be disposed of as electro-scrap.
Consult your municipal authority or your dealer for information about disposal.

**4.3
Specifications**

Name and model	:	MEDISANA FRM Massager
Dimensions (L x W x H)	:	375 x 335 x 110 mm
Weight	:	2 kg
Power supply	:	230 Volts / 50 Hz
Power consumption	:	50 Watt
Article no.	:	88282
EAN Code	:	4015588882821

In accordance with our policy of continual product improvement, we reserve the right to technical and design changes without notice.

5.1 Warranty and repair terms

Please contact your dealer or the service centre in case of a claim under the warranty. If you have to send in the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

1. The warranty period for **MEDISANA** products is three years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

Please do not hesitate to contact us if you have questions about this or any other product.

MEDISANA AG

Itterpark 7-9

D-40724 Hilden

Tel.: +49 (0) 2103 / 2007-60

Fax: +49 (0) 2103 / 2007-626

eMail: info@medisana.de

Internet: www.medisana.de

The service centre address is shown on the last page.

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D&L Marketing Oy
ALERE Electronics AB