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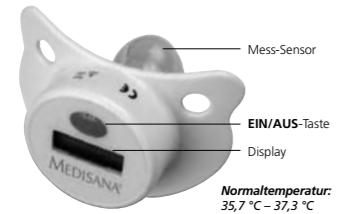
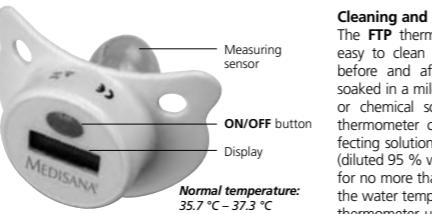
GB) PACIFIER THERMOMETER FTP**Safety instructions**

The **FTP** pacifier thermometer is exclusively intended for measuring the body temperature of children up to the age of 3.

• Under certain circumstances the measurements displayed by the pacifier thermometer may be too low (e.g. if the child is crying or continuously breathing through its mouth). For this reason, in case of doubt or if the child is showing fever-like symptoms (reddening of the skin, sweating, anxiety, refusing to drink), always take a rectal control measurement using a normal fever thermometer if the measurement exceeds 38,5 °C.

• This thermometer is not suitable for continuous use as a pacifier.

• Protect the thermometer from bumps and impacts, and do not drop.

• Avoid ambient temperatures of greater than 60 °C. Do not place the thermometer in boiling water. Do not expose to direct sunlight for long periods.

• The thermometer cannot/should not be opened. The battery cannot be replaced for safety reasons.

• Stop using the thermometer if the teat is damaged.

• Do not attempt to repair the device yourself if problems occur. Stop using the device. Contact your specialist dealer and have all repairs carried out by **MEDISANA** Service.

• Only use normal commercial disinfectant to clean the thermometer, as described in the section entitled „**Cleaning and Disinfecting**“. Do not place damaged thermometers in liquid disinfectant.

• The measuring accuracy does not need to be calibrated during the three-year service life of the battery.

Switching on the thermometer
Switch the thermometer on by pressing the **ON/OFF button** above the display. A short acoustic signal („beep“) indicates that the thermometer has been switched on. A test is performed in which all display segments illuminate for 2 seconds. At an ambient temperature of less than 32 °C an „L“ and a flashing “C“ appear at the top-right-hand side of the display. The thermometer is now ready for operation.

Function test
Every time the device is switched on, the thermometer is checked for correct operation. If a malfunction is detected (inaccurate measurement) “**ERR**” appears in the display and it is no longer possible to take a measurement, the thermometer must be replaced. If you do not switch the device off it will switch off automatically after about 10 minutes.

Use of the thermometer
The pacifier thermometer is only suitable for measuring the body temperature of children by mouth, and must not be used as a pacifier. Put the thermometer in the child's mouth like a normal pacifier. Ensure that the thermometer remains in the child's mouth during measurement, and that the mouth remains closed if possible. The measurement preparation time is 90 seconds. During this time the “C” symbol flashes. When the acoustic signal („beep“) is heard and the “C” symbol stops flashing, the temperature increase has been less than 0,1 °C for 32 seconds. A melody lasting for approx. 10 seconds indicates that the measurement is complete, and the measurement can be read off and stored.

Conditions of guarantee and repair terms
Please contact your dealer or the service centre in case of a claim under the warranty. If you have to send in the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

1. The warranty period for MEDISANA products is three years from date of purchase. In case of a warranty claim, the date of purchase has to be proved by means of the sales receipt or invoice.

2. Defects in material or workmanship will be removed free of charge within the warranty period.

3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.

4. The following is excluded under the warranty:

a. All damage which has arisen due to improper treatment, e.g. nonobservance of the user instructions.

b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.

c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.

d. Accessories which are subject to normal wear and tear.

5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

MEDISANA AG
Iiterrpark 7-9
D-40724 Hilden
Tel.: +49 (0) 2103 / 2007-60
Fax: +49 (0) 2103 / 2007-626
eMail: info@medisana.de
Internet: www.medisana.de

If servicing is required, please contact:

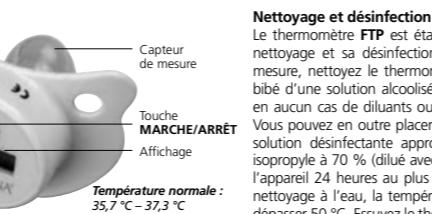
MEDISANA HEALTHCARE UK LTD.
City Business Centre 41 St.
Olai's Court Surrey Quays
London SE16 2XB
Tel.: +44 / 207 - 237 88 99
Fax: +44 / 207 - 252 22 99
eMail: info@medisana.co.uk
Internet: www.medisana.co.uk

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Memory
When the **ON/OFF button** is pressed and the display test has been carried out, the last measurement stored will be automatically displayed for 2 seconds. An “M” for Memory appears in the display at the same time.

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F) THERMOMETRE-TETINE FTP**Cleaning and Disinfecting**

The **FTP** thermometer is watertight and is therefore easy to clean and disinfect. Clean the thermometer before and after each measurement using a cloth soaked in a mild, alcohol-based solution. Do not use thinner or chemical solvents under any circumstances. The thermometer can also be placed in a suitable disinfecting solution. Use a 70% isopropyl alcohol solution (diluted 95% with water). Leave the device in the liquid for no more than 24 hours. When cleaning with water, the water temperature must not exceed 50 °C. Dry the thermometer using a clean cloth. Store the device in a location that is free of dust, dry and hygienic.

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