## Intensiv-Massagegerät ITA

### Intensive massager ITA

Appareil de massage intensif ITA

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- 1 Massage head
- 2 Acupressure attachment
- 3 Point massage attachment
- 4 Power cord with plug
- 5 Slide switch
- 6 - OFF (OFF)
- 7 - Low percussion intensity (LO)
- 8 - High percussion intensity (HI)

*Recommended massage zones:*

- 9 Arms
- 10 Back
- 11 Feet
- 12 Thighs
- 13 Shoulders

It is important to read the safety information carefully before you use the massage device. Keep the instructions for use handy.

Only use the massage device indoors!

Do not use the device over any kind of water container.

### 1.1 Thank you very much

Thank you for your confidence in us and congratulations! By purchasing the **ITA** intensive massage device, you have acquired a high-quality product made by **MEDISANA**. To ensure that you are able to obtain the results you are looking for and can enjoy using your **MEDISANA** massage device **ITA** for a long time to come, we recommend that you read carefully the following information about the use and care of the device.

## 1.2 About the power supply

- Before connecting the device to your power supply, make sure that it is switched off and that the electrical voltage specified on the device's type label matches that of your mains power.
- Keep the power cord and the device itself away from heat, hot surfaces, damp and liquids. Avoid touching the plug or power switch with wet or even damp hands, or when you are standing in water.
- Parts that conduct electricity on the device must not come into contact with liquid.
- If the device falls into water, do not try to pick it up. Instead pull the plug out immediately.
- Always make sure that the plug is easily accessible when the device is connected.
- Disconnect the plug from the socket immediately after using the device.
- Always pull on the plug itself when disconnecting the device from the power supply. Avoid pulling on the power cord!
- Do not carry, pull or rotate the device by its power cord.
- Make sure at all times that the power cord does not represent a tripping hazard. The cord should not be kinked, pinched or twisted.

## 1.3 Specific groups of people

- Particular care should be taken when using the device in the vicinity of children, people with illnesses, and particularly people with special needs.
- Do not leave the device unattended while it is plugged in.
- Do not allow children to use the device as a toy. The device should be stored out of reach of children.
- Do not allow children to play with the clear plastic packaging, as it may pose a risk of suffocation!
- If the device is used from, on or near children or people who are handicapped and can only take limited use of the device, this should always be monitored. Medical devices are no toy!

- Do not use this device if you are pregnant.
- People fitted with a pacemaker, artificial joints or electronic implants should in all cases consult a doctor before using the device.
- You should not use this device if you are affected by one or more of the following illnesses and conditions: circulatory disorders, varicose veins, open wounds, bruises, cracked skin, phlebitis.
- If you feel any pain from use of this device, or find the massage unpleasant, stop the treatment and consult your doctor.

#### **1.4 Before operating the device**

- Check the device and its power cord carefully for any damage before each use. The device should not be used if it is in any way defective.
- Do not use the device if there is any visible damage to the device itself or to the power cord, if the device is not working flawlessly or if it has been dropped or has fallen into water. To avoid any risk of injury send the device to the service point to be repaired.

#### **1.5 Concerning operation**

- The device is intended for private use only and should not be employed for commercial purposes.
- Use the massage device only for the purpose specified in the instructions for use.
- The guarantee is invalidated if the device is misused.
- On no account use any device accessories which are not expressly designed for the purpose by the manufacturer.
- Keep the device away from rain and damp. It should not be operated outdoors.
- It is vital to avoid placing, using or storing the device above containers holding liquids, particularly where the device could fall or be knocked into a bathtub or shower tray.
- Do not use the device in the vicinity of aerosol vapour or at the same time as an oxygen device.

- Do not use the device underneath bedding or cushions. Under no circumstances should you cover the device during use.
- Do not use the device for more than twenty minutes at a time, and allow it to cool down for about thirty minutes before using it again.
- Avoid using the device directly on the head, bones or joints, on the front of the throat, in the genital area or around the kidneys.
- Avoid direct contact between this device and any swollen or inflamed parts of the body or rashes.
- During operation, keep the device well away from the eyes and other sensitive parts of the body.
- It is best not to use this device directly before retiring - the massage that it delivers has a stimulating effect and may delay sleep.
- For understandable reasons, the manufacturer cannot bear the costs of any damage that may result from non-compliance with the instructions for use.
- Should you give this massage device to another person, it is vital that you also pass on these instructions for use.

## 1.6 Maintenance and cleaning

- You the user should not carry out any work on the device, except for cleaning. To avoid possible injury, do not attempt to repair the device yourself. Contact the service point for all repairs.
- Before cleaning the device, switch it off, unplug it from the socket, and allow it to cool. Never immerse the device in water or other liquids.
- Read the information in the section "Miscellaneous, 4.1 Care and storage" on page 19.
- If a defect occurs, do not repair the device yourself. This would automatically invalidate the guarantee. The service point alone should be charged with any repair work.
- When the device is not in use it should be stored in a clean and dry place, ideally in its original packaging.

### 2.1 Scope of delivery and packaging

Please check directly that nothing is missing from the device. The scope of delivery includes:

- 1 **MEDISANA ITA** massage device
- 2 massage attachments
- 1 operating instructions booklet

The packaging can either be reused or recycled. Please dispose of unwanted packaging material responsibly. Should you notice any damage caused in transit when unpacking your device, please contact your dealer promptly.

### 2.2 What is massage?

Massage is one of the oldest methods of treatment known to human kind, and even today is an important element in physical therapy. Mechanical manipulation of the body's outer layers of tissue loosens up tense muscles, boosts circulation and stimulates the regeneration of tissue. Massage can therefore have a stimulatory or relaxing effect, and can furthermore help to alleviate neuralgia and muscular pain.

Massage additionally has a beneficial effect on the vascular system, the internal organs and the body's metabolism.

There are generally held to be five traditional types of massage: Effleurage, kneading, tapotement and vibration massages, as well as rubbing/tapping.

On top of these come such important massage techniques as manual lymph drainage, reflexology, underwater pressure-jet massage, brushing and finger-pressure massage (Shiatsu, acupuncture).

### 2.3 How does the ITA massage device work?

Tapotement is the technique used in this specialised massage device. When set to the lower intensity level, the massage head pulses at a speed of around 2600 pulses per minute, and at 3400 pulses per minute on the higher level, producing a tapotement effect on the area of the body under treatment.

### 3.1 First-time use

The mains plug must be inserted in a suitable socket (230V~) before the device can be used. Please ensure that the device is switched off before plugging it in. The slide switch **5** must be in the *OFF* **6** position.

### 3.2 Operation

Read and note the safety information in Chapter 1 before you start to use the device!

Commence massage by moving the two-step slide switch **5**. The first massage level is activated when you move the slide one step. Set like this, your device emits low intensity vibrations (*LO*) **7**. The higher-intensity (*HI*) massage level **8** is activated when you continue to slide the switch until it catches a second time. To switch the device off, just reverse the switching-on process.

### 3.3 The ideal massage

- You can use this device to massage either yourself or a partner. Switch on the device as described above, and position the massage head **1** on the part of the body you wish to massage. Do not press down, just touch the body lightly with the device. Avoid spending too much time on the same place, rather, you should switch the massage from one area to the other regularly during use.
- You can alter the effect of the massage by altering the direction in which you stroke the device: Massaging towards the heart has a relaxing effect, while massaging away from the heart is energising.
- Massage sessions with the device should last no longer than 10 to 15 minutes. Use the device regularly to achieve the best results.
- After treatment, switch the massage device off, pull out the mains plug and let the device cool down.


### 3.4 Using the attachments

Two massage attachments are supplied with the device - one for acupressure **2** and one for point massage **3**. First switch the device off. Fit whichever attachment you require onto the massage head **1** and restart the device. Carry out the massage as described in section 3.3.

### 4.1 Cleaning and maintenance

- Before cleaning the device, unplug it from the socket, and allow it to cool.
- The device should be cleaned using a soft, slightly dampened cloth, and mild soap if required.
- Do not use harsh cleaning agents or stiff brushes.
- Wait until the device is completely dry before using it again.
- Undo any kinks in the cable.
- You should ideally store the device in its original packaging and keep it in a clean, dry place.

### 4.2 Disposal

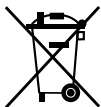This product must not be disposed together with the domestic waste.

All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner.

Consult your municipal authority or your dealer for information about disposal.

## 4.3

**Technical data**

Name/Model : **MEDISANA ITA massage device**

Power supply : 230V~ 50Hz

Power consumption : 25W

Protection category : II 

Item number : 88290

EAN number : 40 15588 88290 6

**We reserve the right to make technical and design changes in the course of continuous product improvement.**

## 5.1 Warranty and repair terms

Please contact your dealer or the service centre in case of a claim under the warranty. If you have to send in the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

1. The warranty period for **MEDISANA** products is three years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
  - a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
  - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
  - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
  - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

### **MEDISANA AG**

Itterpark 7-9  
40724 Hilden

Tel.: +49 (0) 2103 / 2007-60  
Fax: +49 (0) 2103 / 2007-626  
E-mail: [info@medisana.de](mailto:info@medisana.de)  
Website: [www.medisana.de](http://www.medisana.de)

The service centre address is shown on the last page.

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D&L Marketing Oy
ALERE Electronics AB