Medibeam

Akne-Lichttherapie
Acne Light Therapy

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Fototerapia anti-acne

Acne lichttherapie

Valoterapia aknen hoitoon

Ljusterapi mot akne

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**Thank you
very much**

Thank you for your confidence in us and congratulations!
You have acquired a **MEDISANA** quality product with this **Medibeam** acne light therapy unit.

This unit is designed for light therapy of skin diseases, such as acne, pimples, small wounds and herpes simplex.

We recommend that you carefully read the following information on use and care of the product so that you experience its benefits and enjoy using your **MEDISANA Medibeam** acne light therapy unit for many years.

Important

Please read the following safety information carefully before using the device and keep these instructions in a safe place for later use.

**1.1
About the power
supply**

- Before connecting the device to your power supply, make sure that it is switched off and that the electrical voltage specified on the device's type label matches that of your mains power.
- Only connect the mains adapter whilst the unit is switched off.
- Do not operate this unit near high frequency electromagnetic sources.
- Keep the power cord and the device itself away from heat, hot surfaces, damp and liquids. Avoid touching the plug or power switch with wet or even damp hands, or when you are standing in water.
- Parts that conduct electricity on the device must not come into contact with liquid.
- If the device falls into water, do not try to pick it up. Instead pull the plug out immediately.
- Always make sure that the plug is easily accessible when the device is connected.
- Disconnect the plug from the socket immediately after using the device.
- Never pull on the power cord to disconnect the appliance from the mains. Always pull out the mains plug!
- Do not carry, pull or rotate the device by its power cord.

**1.2
Specific groups
of people**

- Particular care should be taken when using the device in the vicinity of children, people with illnesses, and particularly people with special needs.
- Do not leave the device unattended while it is plugged in.
- Children may only use the appliance under adult supervision. Keep the unit out of reach of children. Medical products are not toys!
- The device should not be used by anybody who has any kind of physical disability which could limit their ability to operate the device correctly.

1.3 Before operating the device

- Check the device and its power cord carefully for any damage before each use. The device should not be used if it is in any way defective.
- Do not use the device if there is any visible damage to the device itself or to the power cord, if the device is not working flawlessly or if it has been dropped or has fallen into water. To avoid any risk of injury send the device to the service point to be repaired.

1.4 Concerning operation

- You must wear the protective goggles enclosed with the unit when treating your face.
- Only use the unit as intended and in accordance with the instructions.
- The guarantee is invalidated if the device is misused.
- Use of the **Medibeam** acne light therapy unit is no replacement for a medical diagnosis or treatment.
- Check with your doctor if you have any type of pain or illness before using the unit for the first time.
- If you have any medical questions relating to the use of this unit, your doctor will be able to help you.
- Always seek the advice of a doctor in case of doubt and follow his instructions for treatment.
- If use of the unit is unpleasant or painful, discontinue using it immediately. You may also wish to consult your doctor.
- Persons, whose sensitivity is disturbed or impaired may only apply the device according to the instructions of a physician.
- Keep the device away from rain and damp. It should not be operated outdoors.
- It is vital to avoid placing, using or storing the device above containers holding liquids, particularly where the device could fall or be knocked into a bathtub or shower tray.
- Do not use the device in the vicinity of aerosol vapour or at the same time as an oxygen device.
- Do not use the device underneath bedding or cushions. Under no circumstances should you cover the device during use.
- For understandable reasons, the manufacturer cannot bear the costs of any damage that may result from non-compliance with the instructions for use.
- Should you give this massage device to another person, it is vital that you also pass on these instructions for use.

1.5 Concerning maintenance and cleaning

- You the user should not carry out any work on the device, except for cleaning.
- Before cleaning the unit, switch it off and unplug it from the mains. Never immerse the device in water or other liquids.
- Read the information in the section "Miscellaneous, 4.1 Cleaning and maintenance".
- If a defect occurs, do not repair the device yourself. This would automatically invalidate the guarantee. The service point alone should be charged with any repair work.
- When the device is not in use it should be stored in a clean and dry place, ideally in its original packaging.

2.1 Scope of delivery and packaging

Please check first whether the device is complete.

The package contains:

- 1 **MEDISANA Medibeam** acne light therapy unit
- 1 mains adapter with cable and plug
- 1 spacer ring
- 1 pair of goggles
- 1 instruction manual

The packaging can be reused or recycled. Please dispose properly of any packaging material no longer required. If you notice any transport damage during unpacking, please contact your dealer without delay.

WARNING

To avoid risk of suffocation, keep packaging film away from children.

2.2 How light therapy with Medibeam can help you

Light therapy is a treatment method which uses light beams to treat diseases. This new form of therapy is used with particular success for the convenient and effective treatment of acne. For the **Medibeam** acne light therapy unit, diodes have been developed which are capable of emitting red light of the correct wavelength (about 660 nm). This wavelength promotes and accelerates the healing processes of skin diseases. Therapy with **Medibeam** is safe and painless.

Some of the applications are:

- acne and pimples
- skin infections
- minor lesions and small wounds
- Herpes simplex

2.3 The advantages of light therapy using Medibeam

- The light source of the **Medibeam** acne light therapy unit has special characteristics which accelerate the natural healing processes so that use of the unit may enable the patient to reduce the medication and treatment required.
- The local application of light therapy does not affect any other organs or parts of the body.

3.1 Getting started and using the unit

1. Connect the enclosed mains adapter **6** to the power socket **4** of the unit using plug **5**. Never use a different mains adapter!
2. Before connecting the mains adapter, check that the unit is switched off (switch **3** in position **OFF**).
3. Connect the mains adapter to the correct type of socket (230 V AC). When the unit is ready to operate and functions correctly, the green indicator LED **2** will illuminate. If this LED blinks, then the unit is not working correctly. Do not use the unit in this event and contact the service department.

WARNING

Never look directly into the active light source! You must wear the protective goggles **7 enclosed with the unit when treating your face!**

4. Depending on the intended use, place switch **3** in position **CW** or **PULSE**. **CW** irradiation (CW = constant waves) provides uniform continuous irradiation, while **PULSE** irradiation provides the most intensive light strength and deep pulsed penetration of the beams. In general, treatment of 6 - 8 minutes is recommended for **CW** operation, followed by treatment with **PULSE** operation for the same length of time.
5. To treat the skin with the optimum light intensity, set the distance between skin and light therapy unit so that the smallest possible section of skin is irradiated by the red light. Direct light source **1** precisely onto the part of the skin to be treated.
6. Use of spacer ring **8** to achieve this is recommended. Place spacer ring **8** directly on light source **1** and then place it onto the part of the skin to be treated. Use of the spacer ring guarantees the correct distance between the **Medibeam** unit and the area to be treated and an optimum concentration of the light beams.
7. When you have finished using the unit, place switch **3** into the **OFF** position and disconnect the mains adapter from the mains socket.

NOTES

Hold the light therapy unit vertically above the skin section to be treated.

Avoid direct illumination of the skin area to be treated by other light sources during treatment.

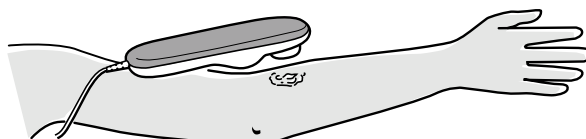
Fig. 1

Small wounds can heal after a few treatments

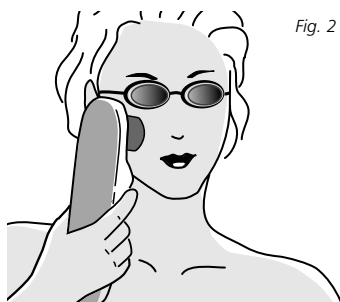
Fig. 2

Face treatment with Medibeam & spacer ring

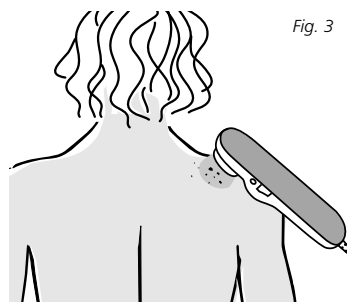
Fig. 3

Treatment of the shoulder using Medibeam

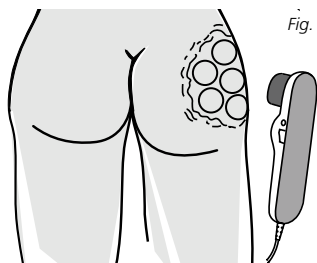
Fig. 4

Larger areas of skin must be divided into smaller sections for treatment

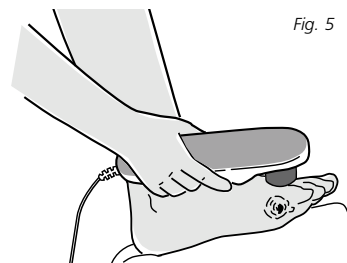
Fig. 5

Support treatment by additional measures (e.g. soft shoes)

Important

WARNING

Irradiation with Medibeam is no substitute for medical treatment. You should consult a doctor before using the unit for the first time. In case of pain or doubt, discontinue treatment and consult your doctor.

3.2 Use for pimples and acne

**) Dr. Earon, Head of Earon Clinics und Dr. Avikam Harel, Pediatric Dermatology Outpatient Clinic, Dana Medical Center, Ichilov Hospital*

In a double blind study* of patients with acne vulgaris, blocked pores, pimples in the face and irritated skin, the following results were obtained after 6 weeks of treatment:

- 20 % were healed completely
- 50 % were healed almost completely
- 20 % displayed considerable improvement
- Only 10 % showed no improvement.

1. Clean the areas to be treated before treatment. Remove any make-up.
2. Irradiate the area to be treated for 6 - 8 minutes with the switch in the **CW** position. The treatment can be repeated up to three times per day at intervals of at least four hours.
3. Follow the instructions of your doctor with regard to medication and treatment in addition to irradiation with the light therapy unit.

3.3 Treatment of small wounds and minor lesions

1. Clean the area to be treated before treatment. Remove any ointments applied and then re-apply after treatment.
2. A small wound which is a few days old can be healed with 6 - 8 treatments.
3. Wounds in areas of the body which are not well supplied with blood require particular attention. In this case preventative measures should be combined with the treatment with the **Medibeam** light treatment unit. The underlying causes must also be removed (e.g. for wounds on a foot: padded shoes). Follow the instructions of your doctor with regard to medication, ointments and treatment in addition to irradiation with the light therapy unit.
4. Irradiate the part of the body to be treated for 6 - 8 minutes with the switch in the **CW** position followed for a further 6 - 8 minutes in the **PULSE** position. The treatment can be repeated up to three times per day at intervals of at least four hours. It is usual to carry out treatment in the morning and the evening. Satisfactory results can also be achieved with one treatment per day.
5. Skin sections of more than 4 - 5 cm diameter must be divided into separate sections. Each section should have a diameter of about 4 - 5 cm (see Fig. 4) and treated as per step 4.

3.4 Use for herpes simplex

Only use continuous radiation (**CW**) when treating herpes simplex. The treatment duration is 8 - 10 minutes, 2 - 3 times per day. Use spacer ring **8** for ease of use with the area to be treated situated in the centre of the ring.


3.5 Switching off the unit

When treatment is completed, switch off the unit by setting switch **CW/OFF/PULSE 3** to **OFF**. Disconnect the mains adapter from the socket.

4.1 Cleaning and maintenance

- Disconnect the mains adapter from the socket before cleaning the unit.
- The device should be cleaned using a soft, slightly dampened cloth, and mild soap if required.
- Do not use harsh cleaning agents or stiff brushes.
- Wait until the device is completely dry before using it again.
- Undo any kinks in the cable.
- You should ideally store the device in its original packaging and keep it in a clean, dry place.

4.2 Disposal

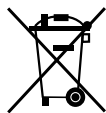Dispose of your old electrical unit in an environmentally acceptable manner when it has reached the end of its service life!

Packaging is re-usable or may be re-cycled; metal parts must be disposed of as scrap metal; plastics, electrical and electronic components must be disposed of as electro-scrap.

Consult your municipal authority or your dealer for information about disposal.

4.3 Directives / Standards

The **Medibeam** acne light treatment unit of **MEDISANA** carries the CE mark (mark of conformity) "CE 0473" which confirms conformity of the unit with "EU Council Directive 93/42/EEC of 14 June 1993 on medical products".

EC importer within the meaning of EC Directive 93/42/EEC is **MEDISANA**.

MEDISANA holds the relevant declaration of conformity.

Equipment classification: Enclosure class II 

4.4 Technical Data

Name and model	:	MEDISANA Medibeam acne light therapy unit
Input voltage	:	230 V AC / 50 Hz
Output voltage	:	24 V AC / 200 mA
Wavelength	:	660 nm
Area of focal point	:	2 cm ²
Continuous radiation (CW)		
Power at focal point	:	18 mW (average)
Maximum power density	:	15 mW/cm ²
Pulsed intensive irradiation (PULSE)		
Peak power at focal point	:	75 mW (average)
Frequency	:	100 Hz
Operative time	:	10 % (average)
Average power	:	7,5 mW
Dimensions (L x W x H)	:	230 x 65 x 49 mm
Weight	:	approx. 0,18 kg (without mains adapter)
Length of mains cable	:	approx. 3 m
Item number	:	47000
EAN number	:	4015588470004

CE 0473

Right reserved to make technical and design changes during the course of continuous product improvement.

5.1 Warranty and repair terms

Please contact your dealer or the service centre in case of a claim under the warranty. If you have to send in the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

1. The warranty period for **MEDISANA** products is three years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

MEDISANA AG

Itterpark 7-9

D-40724 Hilden

Tel.: +49 (0) 2103 / 2007-60

Fax: +49 (0) 2103 / 2007-626

eMail: info@medisana.de

Internet: www.medisana.de

The service centre address is shown on the last page.

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D&L Marketing Oy
ALERE Electronics AB