MEDISANA®

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1.1 Thank you very much

Thank you for your confidence in us and congratulations on your purchase! By purchasing the manicure and pedicure device **Medistyle**, you have become the owner of a high-quality product by **MEDISANA**. Beautiful healthy hands and feet are a matter of good care. Regular treatment – a manicure and pedicure – with the **MEDISANA Medistyle** can improve the quality of your nails, produce soft skin on your hands and feet, and ensure a cultivated and well-groomed appearance. In order to ensure that you achieve the success desired and enjoy using your **MEDISANA Medistyle** for a long time to come, we recommend that you read the following instructions carefully on use and maintenance of the device.

1.2 Important

Please Please read the following safety information carefully before using the device and keep these instructions in a safe place for later use.

- Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your mains supply.
- Only use this unit according to the intended purpose and as directed in the instructions for use. The device is not intended for commercial use.
- Do not use any auxiliary parts that have not been recommended by the maker.
- Do not put any foreign objects into any of the openings. Avoid contact with pointed or sharp objects. The warranty is null and void if it is not used for the intended purpose.
- If the unit is to be used by or on people who are handicapped in any way or in the presence of children, this must always be done under supervision. Medical equipment is not designed to be used as a toy!
- Do not attempt to repair the unit yourself in the event of any defects. Not only does this render the warranty null and void but you can also run the risk of serious hazards (fire, electric shock, injury). Only have repairs carried out by authorised service points.

Please note the following in particular concerning electrical equipment:

- Only connect the mains adapter **1** when the unit is switched off.
- Do not reach for the unit if it has fallen into water. First pull out the mains plug at once.
- Do not use any electrical equipment while in the bath or shower. Do not store or keep any electrical equipment in places where it could fall into the bathtub or washbasin.
- Do not use the unit under the bedclothes or under cushions.
- Keep the mains cable away from hot surfaces.
- Never use the unit in rooms in which aerosols (sprays) are being used or where pure oxygen is being administered.
- Never use the unit if the power cable has been damaged, the mains adapter is defective, if the unit does not work perfectly, if it has fallen into water or has been damaged.
- Never carry, pull or rotate the unit by the mains cable and never pinch the cable between other objects.
- Do not allow the unit to become overheated due to excessively long use.

1.3 Important recommendations concerning your health

- Only pull the plug out of the socket when you have dry hands, when the unit has been switched off and is not being used any longer, and also do this when you have to replace the heads.
- Consult your doctor first before using this unit if you have diabetes or other related illnesses.
- Pregnant women should note the important precautionary measures and the amount of stress they can take, consult your doctor first if necessary.
- Do not treat any parts of the body that have swellings, burns, inflammation, peeling skin, wounds or other sensitive points.
- The treatment should be pleasant. Stop at once and consult your doctor if you feel pain or the treatment seems unpleasant.

2.1 Scope of delivery

Please check first of all that the unit is complete.

The following parts are included as standard:

- **Medistyle** -unit with mains adapter
- 6 accessory heads (Medistyle **S**)
- 11 accessory heads (Medistyle **L**)
- Storage bag
- Manual

If you notice any transport damage during unpacking, please contact your dealer without delay.

2.2 Warning! Special points concerning the MEDISANA Medistyle

- The sapphire heads have an extremely long life and virtually never wear out.
- Uninterrupted treatment should not last longer than 15 minutes. Then switch off the unit for at least half an hour and allow it to cool down.
- Due to the minimal risk of injury, the **Medistyle** is also especially suitable for use by diabetics or older people.
- The **Medistyle** heads should be cleaned with an alcohol solution for hygienic reasons.
- Do not soften the skin beforehand. The **MEDISANA Medistyle** is quick, convenient and can be used everywhere with safety.
- You will find that it is easy to use the **MEDISANA Medistyle** and will be familiar with it after a very short time.

2.3 Frequently asked questions and answers

Why is the Medistyle especially suitable for use with diabetics?

Because the speed is continuously adjustable and the device stops if too much pressure is used, the unit can also be used safely by persons who have reduced sensitivity to pain or are unable to react quickly to pain. Also, the sapphire accessory heads only remove dead skin.

Why are there so many heads?

Each head is designed to solve a particular problem for the hands and feet and therefore especially effective and safe in use.

Can several people use one MEDISANA Medistyle?

Yes. But in any case it is best that each person has its own set of the accessory heads to prevent the transmission of skin diseases.

They can be obtained separately from **MEDISANA AG**.

3.1 Functions

A number of different heads can be used with the **MEDISANA Medistyle**. Switch the device ❹ off ❷. Change the heads by pulling off the head that is no longer required. You can insert the next part by rotating it lightly and overcoming a slight resistance. If you hear a slight click as it clicks into place, then the head is firmly in place. You can also use heads from other nail care units with the aid of the attachment (Medistyle L only) on the **MEDISANA Medistyle**.

Start slowly for each application and increase the continuously adjustable rotating speed ❷ to suit your individual requirements.

This produces a small amount of dust. The sapphire heads only remove dead skin particles with careful grinding. Living skin particles are not removed. Thus there is no risk of injury.

Stop the treatment at once if you find that the removal of hard skin is unpleasant.

3.2 Pedicure

Recommended sequence for a pedicure (care of the feet):

- 1. Coarse-grained sapphire cone**
Art.-No. 85003

To remove hard skin.

- 2. Long sapphire cone**
Art.-No. 85007
(Medistyle L only)

To remove large areas of hard skin and calluses.

- 3. Round sapphire burr**
Art.-No. 85005

To carefully remove corns or warts.

4. Cylindrical burr
Art.-No. 85009

To remove hardened nail surfaces.


5. Coarse-grained sapphire disk
Art.-No. 85010
(Medistyle L only)

To trim thickened nails by filing.

6. Needle burr
Art.-No. 85006
(Medistyle L only)

To carefully remove in-growing nails and to remove small dead portions of the nails.
Only to be used at slow speed!

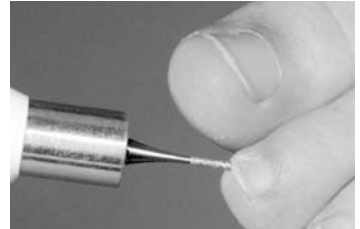3.3 Manicure

Recommended sequence for a manicure (care of the hands):


1. Coarse-grained sapphire cone
(see Pedicure 1)

To remove dry hard skin

2. Flame burr
Art.-No. 85011

To carefully remove in-growing nails.

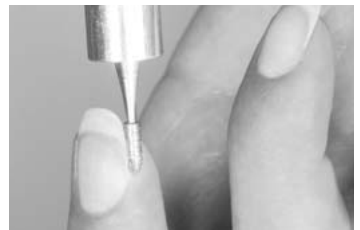3. Fine-grained sapphire disk
Art.-No. 85004


To trim the nails by filing: Hold the unit firmly, run the nail along the disk. Only file in one direction!

4. Device to push back the cuticle
Art.-No. 85008
(Medistyle L only)

To push back the cuticle.

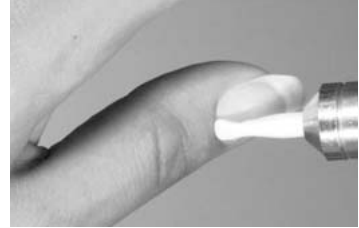
5. Felt cone
Art.-No. 85001

To smooth the edges of the nail and to polish the nail surface.

Attachment
Art.-No. 85015
(Medistyle L only)

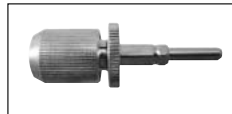To take accessory heads from other units.

Creaming afterwards and regular care with the **MEDISANA Medistyle** contributes to well-cared for hands and feet with soft skin and firm nails with a stable structure.

Any accessory heads that do not come with the **Medistyle S** are available as accessories.

**4.1
Cleaning and
maintenance**

- Always pull the mains adapter out of the socket in the case of an electrical unit.
- Never use strong cleaning agents or stiff brushes.
- Clean the unit with a soft cloth which has been lightly moistened with a mild soap. Never dip the unit into water for cleaning and ensure that no water gets into the unit at any time. Only use the unit again, when it is fully dry.
- Also clean the accessory heads with a small, soft brush if necessary. First allow the heads to dry off completely before you use them again on the unit.
- Store the device in a dry, cool place.
- Fold up the power cable carefully to avoid cable breaks.

**4.2
Disposal**

This product must not be disposed together with the domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

**4.3
Specifications**

Name and model	: MEDISANA Medistyle S / L manicure and pedicure device
Dimensions (Hand-held unit)	: appr. 150 (l) x 40 (Ø) mm
Weight (Hand-held unit)	: appr. 160 g
Current supply	: Operation with mains adapter 230 volts / 50 Hz DC 19 volts, 200 mA
Rotation speed	: 1 600 – 5 000 rpm
Article no. (Medistyle S)	: 85060
Article no. (Medistyle L)	: 85101
EAN-Number (Medistyle S)	: 40 15588 85060 8
EAN-Number (Medistyle L)	: 40 15588 85101 8
LGA/GS-approved safety	

In accordance with our policy of continual product improvement, we reserve the right to technical and design changes without notice.

5.1 Warranty and repair terms

Please contact your dealer or the service centre in case of a claim under the warranty. If you have to send in the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

1. The warranty period for **MEDISANA** products is three years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

MEDISANA AG

Itterpark 7-9
D40724 Hilden

Tel.: +49 (0) 2103 / 2007-60
Fax: +49 (0) 2103 / 2007-626
eMail: info@medisana.de
Internet: www.medisana.de

The service centre address is shown on the last page.

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D&L Marketing Oy
ALERE Electronics AB