Wellnessteine WST

WST wellness stones

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Pietre del benessere WST

Piedras wellness WST

Pedras terapêuticas WST

Wellnesstenen WST

Wellness-kivet WST

Wellness-stenar WST

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1.1 Thank you very much

Thank you for your confidence in us and congratulations on your purchase! By purchasing the **WST** wellness stones, you have gained a high-quality product by **MEDISANA**. To derive the most success and for prolonged enjoyment of your **MEDISANA WST** wellness stones, we recommend that you read the following information on the use and care of the device carefully.

1.2 Important

Please Please read the following safety information carefully before using the device and keep these instructions in a safe place for later use.

- Only use the device for its intended use as in the instructions for use. If misused, all guarantee rights shall become null and void.
- The device is not determined for commercial use.
- Only use the device indoors.
- Do not use any auxiliary parts that have not been recommended by the maker.
- Check the device for any damage before each use (such as cracks, holes etc.). Never use a damaged device.
- Check the stones regularly paying attention to jagged or sharp edges. Only use perfect, rounded stones.
- Be aware that stones might leave scratches on other surfaces (such as wash basins, sinks etc.).
- Do not use any sharp, broken or jagged stones. **There is a danger of injury.**
- When using the tongs to remove stones from the water container, be careful not to drop them.
- Place the device on a firm floor area when in use so that it cannot tip over.
- Never use the device if it is not filled with sufficient water.
- Only fill up the device when the power supply is disconnected.
- Do not fill the water container with hot or boiling water.
- Do not leave the equipment unattended whilst it is connected to the mains.
- Do not put anything in or on the device.
- Disconnect the power supply from the mains after every use and before cleaning the device.
- You may only clean the device yourself. If the device malfunctions do not repair it yourself as this leads to the automatic expiry of all guarantee rights. Only have repairs carried out by authorised service points.
- The device and its attachments are not dishwasher-proof.
- Do not immerse the device in water or other fluids.
- Do not under any circumstances use saucepan scourers, steel wool, strong detergents, paint thinner or alcohol for cleaning the device and its attachments. This can damage the device.

1.2 Important

Please note the following in particular concerning electrical equipment:

- Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your mains supply.
- Do not let the cable hang over the edge of the table, to prevent the device from being accidentally pulled off.
- Keep the device and its power cable away from hot surfaces.
- Do not use the device if it or the power supply cable is damaged, if it does not work properly or if it has been dropped or fallen into water.
- Do not use any electrical equipment while in the bath or shower. Do not store or keep any electrical equipment in places where it could fall into the bathtub or washbasin.
- Do not reach for the unit if it has fallen into water. First pull out the mains plug at once.
- In case of a malfunction, do not try to repair the device yourself. Not only does the warranty expire, but serious danger can result (fire, electric shock, injuries). Only authorised service offices should carry out a repair.
- Do not touch the mains plug when you are standing in water and always take hold of the plug with dry hands.
- Never disconnect the power supply by pulling on the cable. Always pull the power supply unit out of the mains socket.
- Do not use an extension cable.
- Position the mains cable so as to avoid the danger of falling.
- Never carry, pull or rotate the unit by the mains cable and never pinch the cable between other objects.
- Should you give this device to another person, it is vital that you also pass on these instructions for use.

1.3 Important recommendations concerning your health

- If the unit is to be used by or on people who are handicapped in any way or in the presence of children, this must always be done under supervision. Medical equipment is not designed to be used as a toy!
- Do not treat any parts of the body that have swellings, burns, inflammation, peeling skin, wounds or other sensitive points.
- The treatment should be pleasant. Stop at once and consult your doctor if you feel pain or the treatment seems unpleasant.
- If you belong to a risk group due to your health, consult your doctor in all cases before using the wellness stones.
- Do not use the device if you are pregnant or only after consulting your doctor.

2.1 Items supplied and packaging

Please check first of all that the unit is complete.

The following parts are included as standard:

- **WST** water container with insert and lid
- 20 natural stones selected by hand, consisting of:
 - 1 large wellness stone
 - 2 small wellness stones
 - 2 palm stones
 - 6 body stones
 - 8 finger and toe stones
 - 1 massage stone
- 1 pair of tongs
- Manual

The packaging can be reused or recycled. Please dispose properly of any packaging material no longer required. If you notice any transport damage during unpacking, please contact your dealer without delay.

WARNING

To avoid risk of suffocation, keep packaging film away from children!

2.2 How do the MEDISANA WST wellness stones work?

Heat and massage can release tension which is caused by stress, muscle ache and fatigue. The **WST** wellness stones by **MEDISANA** combine massage and heat therapy in a particularly pleasurable way. The use of hot stones is well-known as a healing technique by the Indian and Chinese people. Using massage and heat, certain parts of the body are stimulated which enhances the flow of warmth to the body's energy centres.

The body is flooded with a healthy warmth. These soothing heat rays relieve stress, alleviate tension and loosen up muscles which prepare for a follow-up massage.

For this purpose the **WST** wellness stones are heated up in water and then placed on the energy centres of the body, such as the forehead, stomach, back, hands, feet etc.

The hot stones have many different sizes and shapes and can be used in various ways. You can also use hot and cold stones alternately. This stimulating change of temperature also has a positive effect on the energy channels (Meridians) in your body.

Make **WST** wellness stones by **MEDISANA** an essential part of your wellness routine. Take your time. Prepare your body with a bath or meditation. Choose a quiet, warm and pleasant room. It's even more enjoyable with a friend - now is the time to spoil yourself.

3.1 Getting started

1. Place the device on a level surface that is resistant to water. As a precaution, put a protective sheet underneath to prevent any damage to the delicate surface.
2. Remove the stones from the packaging and wash them one by one. Pay particular attention to jagged edges, sharp or broken stones. **There is a danger of injury.** These stones may also scratch wash basins or sinks.

CAUTION

Make sure you clean the water container and the stones before each use, as well as changing the water.

3. First place the insert in the container to protect the water container. Then carefully place the stones in the water container one by one **2**.
4. Fill up the container with water until all the stones are covered. The water must not exceed or come below the mark for the maximum and minimum level **7**. Do not put in hot or boiling water. This will damage the device. Put the lid **1** on the water container.
5. Plug the power supply into the mains socket. The green **4** and red **5** LED indicators light up. Turn the dial **3** to the required position.
6. Once the set temperature on the dial is reached, the red LED light goes out. It takes approx. 45 minutes to reach the maximum temperature setting for water at 20 °C.
7. The stones have reached the required temperature. Using the tongs supplied **6** take the stones carefully out of the container. Be careful not to drop the stones when using the tongs as they could cause injury.

3.2 General information


1. Before you use the wellness stones, check their temperature with your fingers. If the stones are not too hot, put them in your hand to check they are pleasantly warm in your hand. Leave the stones to cool down briefly.
2. If you do not want to use the stones when wet, dry them with a clean cloth.
3. Put stones that have cooled off back in the water container so they can be warmed up again within a few minutes.
4. When you have finished, disconnect the power supply from the mains socket.
5. Take out all the stones one by one and dry them with a cloth.
6. Empty the water container and dry it carefully.

3.3 Arranging the individual stones

Once you have prepared the stones you can distribute them over your body.

You can do this in any order you wish and for as long as you wish. Proceed carefully at first and try out different combinations.

Cold stones also have a relaxing effect. First place a warm stone on a Meridian point and then a cold one in the same place. This can have a really positive effect on the energy flow in the body.

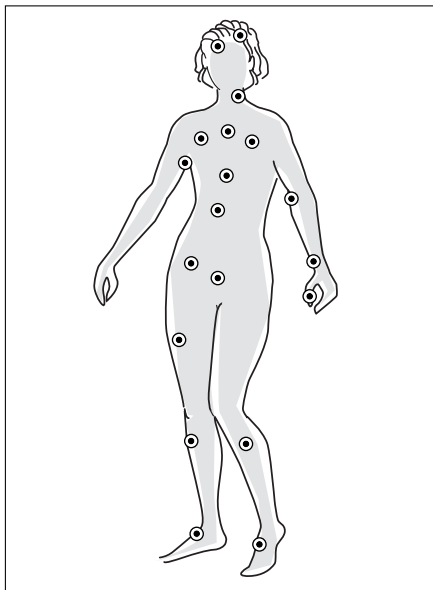*Effective points for
the best arrangement
of the stones.*

Finger and toe stones

These eight stones are intended for placing on fingers and toes. Spread the fingers or the toes. Wedge the stones between the individual fingers or toes. Lie down comfortably on your back to do this. This application can also be used well in combination with others. For instance, place one stone between both foot bones of the large and the second toe and massage this point.


Body stones

These six stones are intended for using on certain parts of the body, such as the neck, shoulders, spine, face etc. They can help to alleviate tension in these areas.

Examples:

Lie down on your stomach, relax and place the warm body stones along the spine. Leave them there until they become cold. This is ideal as preparation for a massage, then a friend can remove the stones and continue your wellness routine with a massage. However, you could always replace the warm stones on the back or at the same time have someone massage you with the large stone in certain points.

Lie down on your back, relax and place the body stones on the upper part of your body. Place the stones in a row. Close your eyes, feel the energy from the stone and enjoy the warming power of nature. This application has a balancing and calming effect and relieves tension.

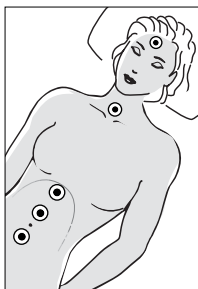Palm stones

These two stones are intended for using in the palms. They then transfer the natural power of the stone to your palms. This works best in the lotus position. Sit comfortably with a straight back. Do not lean against anything. Sitting with crossed legs is ideal for meditating at the same time. Place your hands on your knees or thighs. Your palms are facing upwards. Try to sit still in this position for up to 10 minutes for the first application. With practice you will be able to meditate as long as you want. The effect of this application can be increased by relaxation exercises, such as inhaling and exhaling deeply and/or tensing and relaxing certain parts of the body.


Massage stone

This stone is intended for massaging the body. With a little massage oil the massage stone glides smoothly over the body. First test the oil on a small area of skin for any allergies or irritations before using on the whole body. More heat can be retained for a long time due to the size of the stone. Therefore, it is ideal for massaging larger regions of the body.

For instance, after a tiring day at work or playing sport, massage your arms and legs with the large massage stone in gentle up and down movements.

The massage stone is particularly suitable for treating cellulite on the thighs. Massage the front and back thigh alternately with the stone using gentle up and down movements. This type of massage can also be used to treat cellulite on other parts of the body.

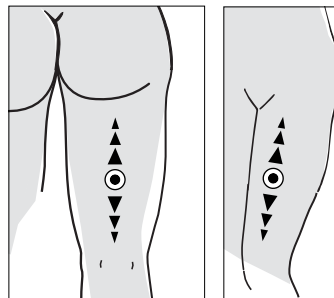Two other particular points on the body which are especially responsive to the treatment of cellulite are shown on the right.

Fig. above: on the bottom:

Place a heated stone on both sides of the base of the spine directly above the middle of the bottom. Massage these points individually if possible with the large massage stone. Get a friend to help if you need it.


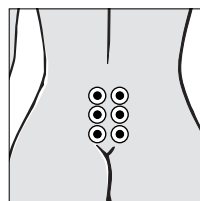
Fig. below: in the middle of the buttock:

The massage stone can be used here as described above.

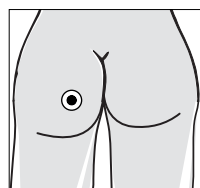
Large wellness stone

The large wellness stone is intended for treating specific parts of the body. Its size enables it to retain the heat for a long time and is especially suitable for large parts of the body. Place it or keep it on the chosen part of the body.

The large wellness stone is like the massage stone and is especially suitable for treating cellulite. Heat the individual problem zone with the large wellness stone and massage gently. Then massage the whole area with the massage stone.

The large wellness stone can be used on any part of the body and wherever you feel it is comfortable. One particular point which is responsible for the the body's sensitivity to cold is located on the inside of the lower leg underneath the knee joint. Massage this point gently with the large wellness stone. A comforting feeling of warmth flows through your body after successfully massaging this point.

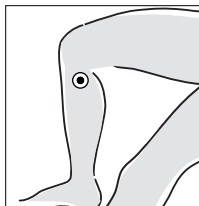
Small wellness stone


Both small wellness stones are intended for using on the face, stomach, back and legs. Their small size make them especially good for using on smaller and uneven parts of the body.

The two small wellness stones can be used for the face with other body stones if required.

One particularly sensitive point is located between the eyebrows and promotes a general feeling of well-being if stimulated by the warm stones. Lie down on your back and let the warm stones radiate inside you.

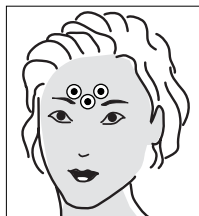Both small wellness stones can be used to help prevent wrinkles from forming on the forehead. Place the stones on both sides of the forehead, one finger's breadth over the eyebrow.

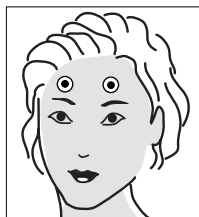4.1 Cleaning and maintenance

- Disconnect the power supply unit from the mains before you clean the device and empty the water container.
- Clean the device and the stones using a soft damp cloth and a mild detergent after every use. Never immerse the device in water to clean it. Renew the water every time you use it.
- Never use abrasive detergents or strong brushes.
- Dry the device and stones properly.
- Store the device and stones in a dry, cool place.

4.2 Disposal

This product must not be disposed together with the domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner.

Consult your municipal authority or your dealer for information about disposal.

4.3 Specifications

Name and model	:	MEDISANA WST Wellness stones
Current supply	:	230 - 240 V~ 50 Hz
Max. power consumption	:	approx. 83 W
Dimensions L x W x H	:	approx. 310 x 210 x 185 mm
Weight	:	approx. 2,3 kg with stones
Item number	:	88150
EAN number	:	40 15588 88150 3

In accordance with our policy of continual product improvement, we reserve the right to technical and design changes without notice.

5.1 Warranty and repair terms

Please contact your dealer or the service centre in case of a claim under the warranty. If you have to send in the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

1. The warranty period for **MEDISANA** products is three years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

MEDISANA AG

Itterpark 7-9
D40724 Hilden

Tel.: +49 (0) 2103 / 2007-60
Fax: +49 (0) 2103 / 2007-626
eMail: info@medisana.de
Internet: www.medisana.de

The service centre address is shown on the last page.

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D&L Marketing Oy
ALERE Electronics AB