EN Explanation of symbols

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
Important Information - Retain for future use!

EN

- **Please read this manual carefully before using this device.**
- **This device is not intended for medical use in hospitals. For household use only!**
- **To avoid any risk of electric shock when using the pad please observe the following:**
 - Do not use if wet. If the pad has become damp or wet, completely dry the pad before use.
 - Only use the pad in dry rooms, not in a bathroom etc.!
 - Check that the pad is undamaged each time before use, i.e. do not use it if there is visible damage to the hand control, switches, plug or the electric heating pad itself.
 - Do not insert pins or other sharp objects.
- **To avoid overheating:**
 - Do not leave the pad switched on unsupervised; switch it off when you leave the room.
 - Do not use the pad folded.
 - Avoid creasing the pad.
 - Do not cover the hand control. To ensure the hand control can function correctly, do not place the hand control on or under the flexible section of the appliance, do not tuck the hand control under pillows, cushions or clothing.
 - The cable should not be wound around the hand control.
- **Do not use on a person incapable of movement, a baby a person insensitive to heat or animals.**
- **Excessive application at a high setting can result in skin burns.**
- **When the appliance is likely to be used for an extended period, for instance by the user falling asleep, controls should be adjusted to a setting recommended for continuous use.**
- **Do not use on a baby. Do not let children play with the electric heating pad.**
Do not leave children unsupervised with the electric heat pad because they cannot correctly recognise or appreciate possible dangers.
- **Examine the device frequently for signs of wear or damage.** If there are signs of wear or damage or the device has been misused, please contact the manufacturer or authorised agent (see below) prior to any further use.
- **If the supply cord is damaged, it must be replaced by the manufacturer, an authorised agent or similarly qualified persons in order to avoid any hazards.**
- **Do not attempt to open or repair the device your self. Special tools are required to re-assemble the device.**
- **When not in use, allow the pad to cool down before folding. Store the electric heating pad flat in a dry place.** Do not place heavy objects on it and never store the pad with sharp creases. Otherwise the pad may be damaged.

1. How to Use a Microlife Heating Pad

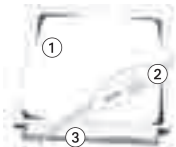- ① Soft fabric cover – removable and machine washable. Press-stud fastener.
- ② Hand control with slider for heat settings and ON indicator.
- ③ Pad

! WARNING: Never use the heating pad without the cloth cover.
Allow some air circulation between the heating pad and skin during use.

1. Put the cloth cover ① on a flat surface. Open the press-stud.
2. Insert the pad ③ completely into the cover, fitting one side at a time. Make sure the cable and hand control ② remained outside the cover.

3. The pad will lie flat, filling the entire cover. Close press-stud.
4. Plug into a 230-240 volt AC mains socket and set hand control to the desired heat setting using the slider.
 - Select the temperature setting on the slider – this switches on the heat pad. An indicator will illuminate, showing that the appliance is connected to the power supply and is switched on.
 - Setting 1 is recommended for continuous use.
 - The electric heating pad is switched off by moving the slider to the «0» position, the indicator will then go out.
5. The electric heating pad hand control will automatically turn off the heater after approximately 90 minutes. However, the indicator lamp will remain illuminated. To reset the timer and start a further heating period, move the slider to the off «0» position and then back to the desired heat setting.
6. Unplug when finished.

2. Cleaning Instructions

! WARNING: Always unplug the power cable from the mains socket before you undertake any cleaning of the product.

1. The detachable cover may be machine washed after removal from the pad according to the washing instructions given on the cover.
2. The pad must not be machine washed or dry-cleaned!
3. Do not allow the hand control to get wet. Do not use any detergents or solvents, they may damage the device.
4. On no account hold the hand control under water. Otherwise moisture may enter the hand control and cause damage.
5. Check the hand control for signs of damage, dry the pad and its fabric cover completely before switching on again.

3. Guarantee

This product is guaranteed for **2 years** from the date of purchase. This guarantee covers the electronic hand control, cable, and the pad itself. This guarantee does not cover incorrect use, accidents or damage resulting from failure to follow instructions, or modifications to the device not carried out by the manufacturer.

The guarantee is only valid upon presentation of the guarantee card filled out by the dealer.

4. www.microlife.com

Detailed user information about our products and about our services can be found at www.microlife.com.

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