## **DIGITAL BABY SCALE 2010**

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# 1 Before initial use

## 1.1 Intended Purpose

The Digital Baby Scale 2010 is the ultimate scale for weighing safely your newborn. This easy to use, large LCD display scale is battery-operated scale and makes daily weighing very easy.

## 1.2 Safety advice

- Always read the safety instructions carefully.
- Follow the basic safety precautions for all electronic equipment when using this Scale.
- Keep this User's Manual for future reference.
- This unit is not waterproof. Don't expose this product to rain or moisture.
- Lay this equipment on a reliable flat surface before setting it up.
- If any of the following situation arises, get the equipment checked by a service technician:
  - The equipment has been exposed to moisture.
  - The equipment has been dropped and damaged.
  - The equipment shows signs of breakage.
  - Do not place the basic unit in a damp room at a distance of less than 1.5 m away from a water source.
- Don't disassemble this scale. The scale doesn't contain any parts that can be repaired by the user.
- Use the scale only for its intended purpose.

## 1.3 Technical specifications

Operation temperature: 10°C ~ 40°C

Environmental humidity: <80% R.H.

Max. capacity: 20kg

Min. capacity: 20g


Graduation: 5g

Power supply: 4 x Non - rechargeable AA batteries (not included)

Unit change: kg or lb

# 2 Installation


Put batteries in the battery compartment to activate the display :

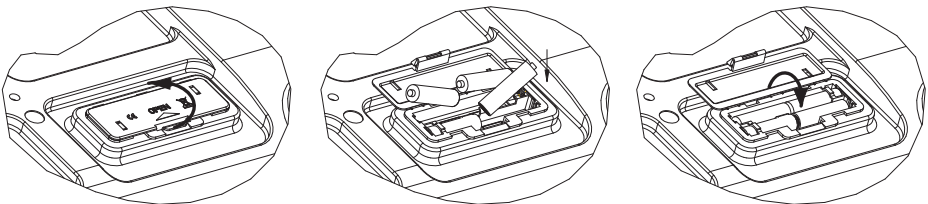- Open the battery compartment at the bottom of the unit.
- Insert the 4 enclosed AA non - rechargeable batteries (1,5V) following the polarity as shown.
- Close the battery compartment.

**Do not short-circuit & dispose of in fire. Remove the batteries if this device is not going to be used for a long period.**

### 3 Description

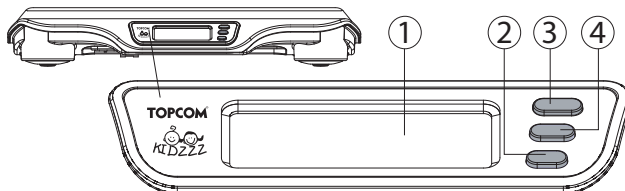1. LCD display
2. MEMORY - button
3. ON/OFF - TARA Button
4. MODE - Button

### 4 Simple operation

- Press the **ON/OFF - TARA** - button (3). The LCD display will show all segments for 2 seconds before entering weight mode.
- Press the **MODE** - button (4) to select 'Kilogram (**kg**)' or 'Pound (**lb**)'
- Gently place the baby on the pan.
- The babies weight (maximum 20Kg) flashes on the display after a few seconds.
- If the baby remain lying on the pan after its weight is shown, the value will continue to flash for 5 minutes and then turn off automatically.

**If you use a towel to lay the baby on, it should be placed on the pan before switching on. If you forgot to do so, press the ON/OFF - TARA button to eliminate the weight of the towel immediately after you have layed down the towel.  
If you remove the towel, the scale will display a negative reading.**

- Press and hold the **ON/OFF - TARA** - button (3) for a few seconds to turn the scale Off or it will turn off automatically in idle mode after 60 seconds.

### 5 Memory

You can store the last result in the memory for later review.

- When the weight flashes on the display, press the **MEMORY** - button (2). '**SET**' will be briefly displayed.

To recall the weight:

- While the scale is in "**0.00**" situation, press the **MEMORY** - button (2).
- The stored weighing value saved in the memory and the '**MEM**' icon is displayed.
- Press the **MEMORY** - button (2) again to leave the memory function.

### 6 Error messages

<b>Lo</b>	Low power - Replace the batteries
<b>Err</b>	Error during initialisation or measurement. Turn the scale off and on again and make a new measurement.
<b>O_Id</b>	Overload. The weight of the baby is over 20kg.

## 7 Cleaning

Clean the scale with a slightly damp cloth or with an anti-static cloth. Never use cleaning agents or abrasive solvents.

## 8 Disposal of the device (environment)

At the end of the product lifecycle, you should not throw this product into the normal household garbage but bring the product to a collection point for the recycling of electrical and electronic equipments. The symbol on the product, user guide and/or box indicate this. Some of the product materials can be re-used if you bring them to a recycling point. By re-using some parts or raw materials from used products you make an important contribution to the protection of the environment. Please contact your local authorities in case you need more information on the collection points in your area.

**Batteries must be removed before disposing of the device.**

**Dispose of the batteries in an environmental manner according to your country regulations.**

## 9 Topcom warranty

### 9.1 Warranty period

The Topcom units have a 24-month warranty period. The warranty period starts on the day the new unit is purchased.

Consumables or defects causing a negligible effect on operation or value of the equipment are not covered.

The warranty has to be proven by presentation of the original purchase receipt, on which the date of purchase and the unit-model are indicated.

### 9.2 Warranty handling

A faulty unit needs to be returned to a Topcom service centre including a valid purchase note.

If the unit develops a fault during the warranty period, Topcom or its officially appointed service centre will repair any defects caused by material or manufacturing faults free of charge.

Topcom will at its discretion fulfill its warranty obligations by either repairing or exchanging the faulty units or parts of the faulty units. In case of replacement, colour and model can be different from the original purchased unit.

The initial purchase date shall determine the start of the warranty period. The warranty period is not extended if the unit is exchanged or repaired by Topcom or its appointed service centres.

### 9.3 Warranty exclusions

Damage or defects caused by incorrect treatment or operation and damage resulting from use of non-original parts or accessories not recommended by Topcom are not covered by the warranty.

The warranty does not cover damage caused by outside factors, such as lightning, water and fire, nor any damage caused during transportation.

No warranty can be claimed if the serial number on the units has been changed, removed or rendered illegible.

This product is in compliance with the essential requirements and other relevant provisions of the directive 89/336/EEC.

The Declaration of conformity can be found on:

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