
waterpik®
treat yourself better®

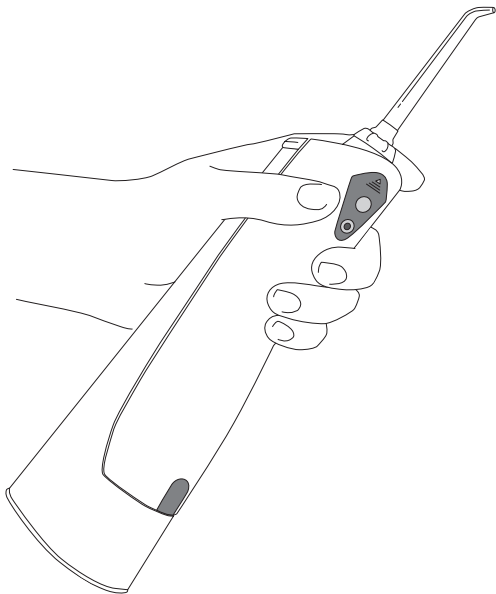Waterpik® Cordless Dental Water Jet

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Irrigador bucal inalámbrico Waterpik®

.....
www.waterpik.com

Congratulations! By selecting the Waterpik® Cordless Dental Water Jet you are now on your way to better gum health, and a healthier smile!


With a clinically proven combination of water pressure and pulsations, the Waterpik® Dental Water Jet removes harmful bacteria deep between teeth and below the gumline – where brushing alone won't reach.

Thank you for choosing the Waterpik® Dental Water Jet.

Traditional Oral Care


Traditional brushing can't easily reach hidden food debris and harmful bacteria.

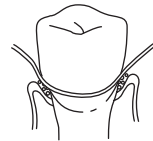Flossing has limited access below the gumline.

Oral Care With The Waterpik® Dental Water Jet

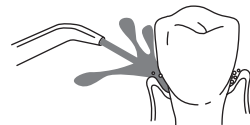The Waterpik® Dental Water Jet removes trapped food debris and harmful bacteria.

Important Safeguards	4
Product Description	6
Getting Started	8
Care of Your Waterpik® Cordless Dental Water Jet	13
Troubleshooting Guide	14
Limited Two-Year Warranty	15

IMPORTANT SAFEGUARDS

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, the basic safety precautions below should always be followed.

READ ALL INSTRUCTIONS BEFORE USING.

DANGER:

To reduce the risk of electrocution:

- Do not handle charger with wet hands.
- Do not immerse in water or other liquid.
- Check the charger cord for damage before the first use and during the life of the product.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub, shower stall or sink.
- Do not reach for a product that has fallen into water. Unplug immediately.



WARNING:

To reduce the risk of burns, electric shock, fire, serious injury or damage to tissue:


- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. Contact Water Pik, Inc.
- Closely instruct and supervise children and individuals with special needs in the proper use of this product.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas. See operating instructions for correct usage.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Only use jet tips that are recommended by Water Pik, Inc.
- Fill reservoir with water or other dental professional recommended solutions only (i.e. mouthwash or prescribed solutions).

- Do not drop or insert any foreign object into any opening on the product.
- Keep charging unit and cord away from heated surfaces.
- Do not operate where oxygen or aerosol sprays are being used.
- Do not use mouthwashes/rinses that contain IODINE or Bleach.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.


SAVE THESE INSTRUCTIONS.

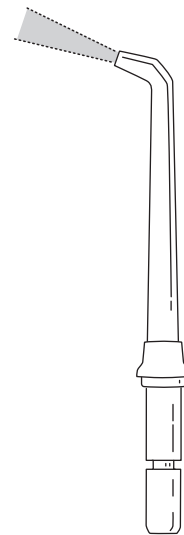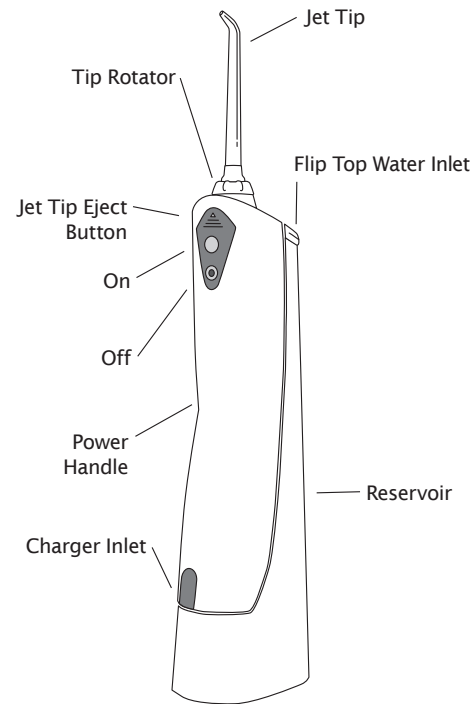Deep Cleaning

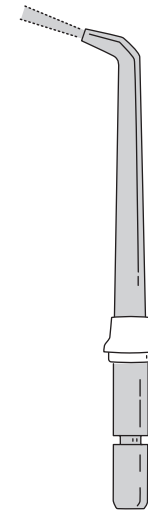Gentle Cleaning

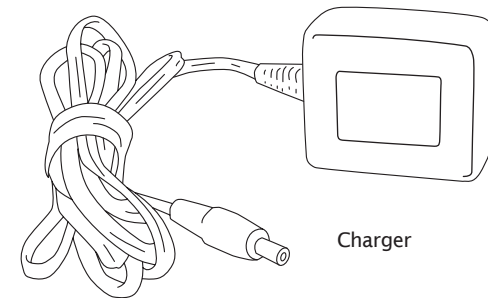
Charger

GETTING STARTED


Important! Turn Dental Water Jet off and charge for 24 hours before first use.

Plug the cord into any standard electrical outlet. Make sure you push the switch into the 'OFF' position before charging. The unit will not charge if the switch is in the 'ON' position. See "Care of your Waterpik® Cordless Dental Water Jet" for guidance on how to look after your battery. The Waterpik® Cordless Dental Water Jet will stay charged for up to 2 weeks depending on use.

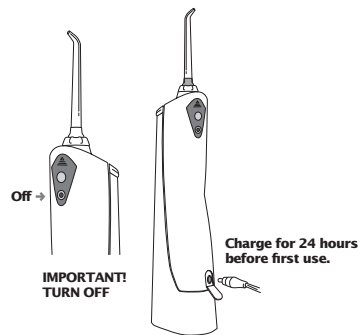
Preparing the Reservoir

Slide Dental Water Jet power handle into back of reservoir and press firmly to snap into place.


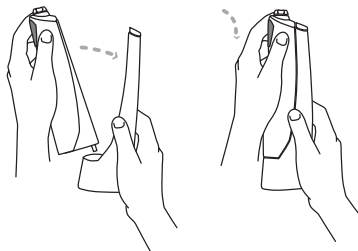
Lift flip top on reservoir and fill the reservoir with lukewarm water. Disconnect charging cord from the unit prior to filling the reservoir and product use.

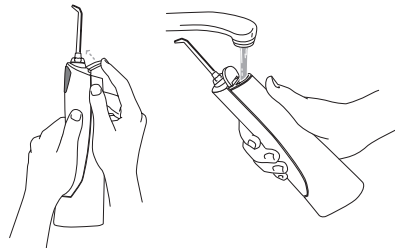Inserting and Removing Tips


Insert tip into the center of the knob at the top of the Dental Water Jet handle. The colored ring will be flush with the end of the knob if the tip is correctly locked into place.

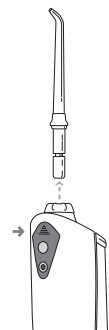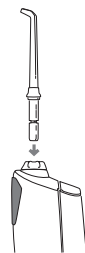To remove the tip from the handle, press the jet tip eject button and pull tip from the handle.

Adjusting the Pressure Setting

You can adjust the pressure setting on your Cordless Dental Water Jet by simply changing jet tips. For low pressure utilize the dark blue tip. For high pressure utilize the clear tip.

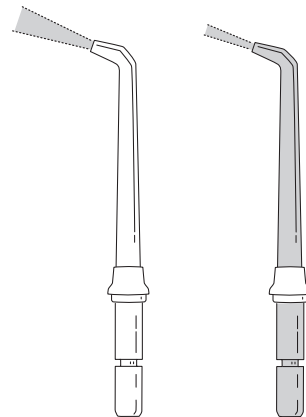
Deep Cleaning


Gentle Cleaning

Placement in Mouth


Lean low over sink and place the jet tip in mouth. Aim the tip toward teeth, and turn the unit on.

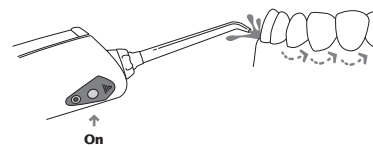Keep unit upright during use for best results.


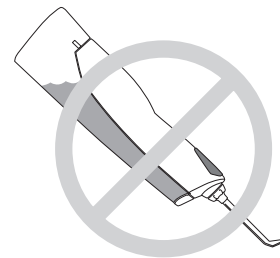
Direct the stream by turning the knob at the top of the handle.

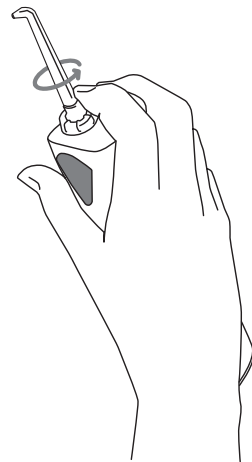
GETTING STARTED

Recommended Technique


Direct the jet stream at a 90-degree angle to your gumline. Close lips to avoid splashing but allow water to flow freely from mouth into the sink. For best results, start in the molar area (back teeth) working toward the front teeth. Glide jet tip along gumline and pause briefly between teeth. Continue until all areas around and between the teeth have been cleaned.

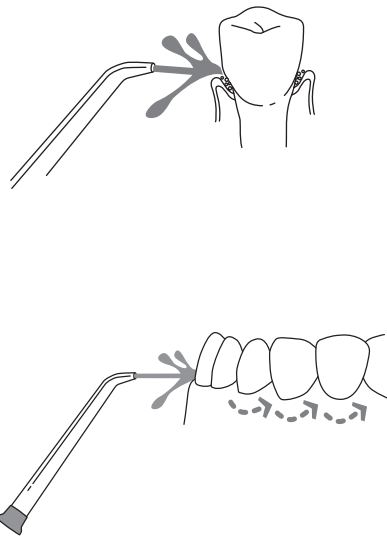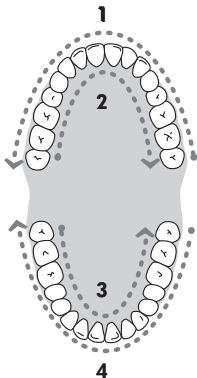**When Finished**

Turn the unit off. Empty any liquid left in the reservoir.

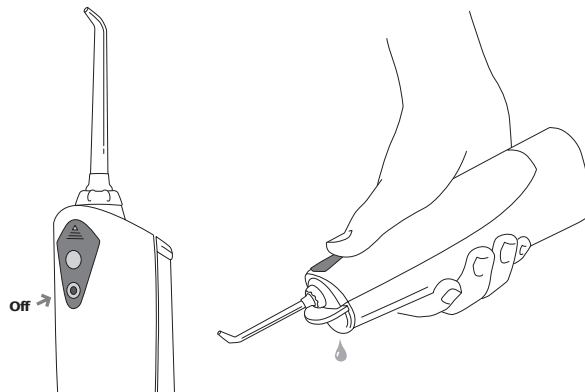If desired, the reservoir may be removed from the power handle by pushing down on the top of power handle and pulling it away from the reservoir.

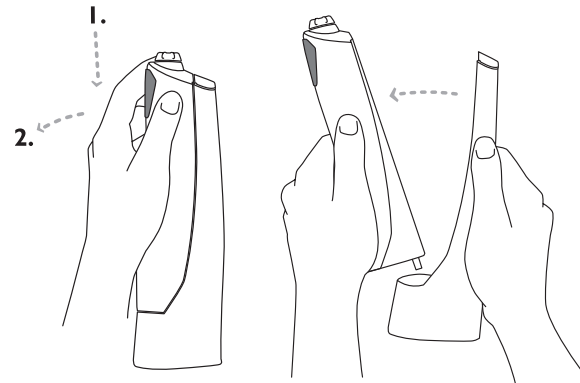GETTING STARTED

Using Mouthwash and Other Solutions

Your Waterpik® Dental Water Jet can be used to deliver mouthwash and antibacterial solutions. Ask your dental professional to recommend the appropriate solution for you. Note: use of some antibacterial solutions could shorten the life of your Waterpik® Dental Water Jet.

After using any special solution in the reservoir, take a moment after every use to rinse the unit. Partially fill the reservoir with warm tap water. Aim the jet tip into your sink. Turn on the unit and let it run until empty. This will help keep the tip from clogging and help your unit to last longer.

CARE OF YOUR WATERPIK® CORDLESS DENTAL WATER JET

Never immerse the unit in water. Wipe the exterior when necessary with a soft cloth and mild non-abrasive cleanser. Reservoir may be washed with soapy water or in top rack of dishwasher.

Before cleaning, unplug from the electrical outlet.

Before exposing the Dental Water Jet to freezing temperatures, remove the reservoir and run the unit until completely empty.

Battery Care

Before first use, charge for 24 hours to obtain full charge.

To charge your Cordless Dental Water Jet, plug the cord into any standard electrical outlet. Make sure you push the switch into the 'OFF' position before charging. The unit will not charge if the switch is in the 'ON' position.

If you use your Cordless Dental Water Jet once a day or less, it is recommended you charge it overnight, once a week. If you use it twice a day or more, charge the unit more frequently. Do not run the batteries down completely as this may shorten the lifespan of the batteries and will require charging the unit for 24 hours to obtain a full charge.

If the Cordless Dental Water Jet is to be stored for an extended period of time, charge it for 24 hours before storing. Never store the Cordless Dental Water Jet for extended periods unless it is fully charged. To maximize the lifespan of the batteries the battery charger may be kept connected to the unit continuously, except for when the product is in use.

Battery disposal

The unit contains a Nickel Metal Hydride battery. At the end of the battery life, please recycle the unit at your nearest recycle center per local or state requirements. The battery is not replaceable.

Service Maintenance

Waterpik® Dental Water Jets have no consumer-serviceable electrical items and do not require routine service maintenance. For all your service needs go to www.waterpik.com. Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit.

When outside the United States or Canada: Contact your local dealer or write Water Pik, Inc. for the address of the nearest authorized service center.

TROUBLESHOOTING GUIDE

PROBLEM	CAUSE	CORRECTIVE ACTION
Jet Tips Leakage between the tip and the tip rotator.	1) Tip not fully engaged.	1) Remove tip and reinstall.
Dental System Reservoir leaks.	1) Liquid was left in reservoir. 2) Power handle not correctly engaged with reservoir. 3) Water inlet is open.	1) Always empty the reservoir after use. 2) Remove reservoir then reinstall as in instruction steps 6 and 1. 3) Close water inlet.
Inadequate pressure/No fluid flow.	1) Unit is upside down. 2) Insufficient liquid in reservoir. 3) Reservoir connection clogged.	1) Jet tip must be higher than reservoir for correct water flow. 2) Fill with liquid. 3) Fill reservoir with warm tap water and one cap of vinegar. Run until empty. Rinse with full reservoir of tap water.
Unit doesn't start.	1) Battery not charged.	1) Be sure outlet is functioning; test with another appliance. See "Care of your Waterpik® Cordless Dental Water Jet" for details on how to charge the battery.

Still have questions? Go to www.waterpik.com.

LIMITED TWO-YEAR WARRANTY

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 2 years from the date of purchase. Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will replace at its discretion any part of the product, which in its opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. Installation is the responsibility of the consumer and is not covered by the warranty. This limited warranty excludes accessories or consumable parts such as jet tips, etc.

This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.

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Water Pik, Inc.

1730 East Prospect Road
Fort Collins, CO 80553-0001 USA

www.waterpik.com

Waterpik® (stylized) is a trademark of Water Pik, Inc. registered in the European Community, Switzerland, Mexico, Canada, Japan, and the United States.

Waterpik® is a trademark of Water Pik, Inc. registered in Argentina, India, Spain, Korea, Canada, Brazil, United Kingdom, Ukraine, Chile, France, Finland, Benelux, Australia, Switzerland, Taiwan, South Africa, Sweden, Norway, Russian Federation, Czech Republic, New Zealand, China, Mexico, Austria, Hungary, Singapore, Italy, Germany, Denmark, Poland, and the United States.

Treat Yourself Better® is a trademark of Water Pik, Inc. registered in the United States.

Designed in USA.
Printed and assembled in CHINA.

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Treat yourself better® ist eine in den Vereinigten Staaten eingetragene Marke der Water Pik, Inc.
Treat yourself better® est une marque de commerce de Water Pik, Inc. déposée aux États-Unis.
Treat yourself better® è un marchio registrato di Water Pik, Inc. registrato negli Stati Uniti.
Treat yourself better® is een handelsmerk van Water Pik, Inc. geregistreerd in de Verenigde Staten.
Treat yourself better® es una marca de comercio de Water Pik, Inc., registrada en Estados Unidos.